

# STEMprep™ Mouse Spleen Dissociation Kit



**For processing 50 spleens using the STEMprep™ Tissue Dissociation System**

Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

Catalog #100-2138

50 Preparations

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

STEMprep™ Mouse Spleen Dissociation Kit provides an optimized enzymatic dissociation cocktail for processing spleen tissues into single-cell suspensions. When used with the STEMprep™ Tissue Dissociator (Catalog #100-2112), mechanical dissociation is combined with enzymatic degradation of the extracellular matrix, efficiently breaking down tissue structure while preserving cellular integrity. This approach yields a high number of leukocytes, particularly dendritic cells, while preserving cell surface epitopes. The resulting single cells are immediately suitable for downstream applications, such as cell separation, culture, or various analyses.

For best results, use with STEMprep™ Tissue Dissociator and STEMprep™ Sample Tubes (Catalog #200-0800).

## Product Information

The following components are sold as part of STEMprep™ Mouse Spleen Dissociation Kit (Catalog #100-2138) and are not available for individual sale.

COMPONENT NAME	COMPONENT #	SIZE	STORAGE	SHELF LIFE
STEMprep™ Enzyme A	100-2131	1 x 2.5 mL	Store at -20°C.	Stable until expiry date (EXP) on label.
STEMprep™ Enzyme E	100-2160	1 x 1.25 mL	Store at -20°C.	Stable until expiry date (EXP) on label.
STEMprep™ Enzyme Diluent Z	100-2134	1 x 125 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.

## Materials Required but Not Included

PRODUCT NAME	CATALOG #
50 mL Conical Tubes	e.g. 38010
70 µm Cell Strainer	e.g. 27260
Culture Dish, Non-Treated	e.g. 38045
D-PBS (Without Ca++ and Mg++)	37350
EasySep™ Buffer OR PBS containing 2% FBS and 1 mM EDTA	20144
STEMprep™ Sample Tubes	200-0800
STEMprep™ Tissue Dissociator	100-2112

## Preparation of Reagents and Materials

Thaw STEMprep™ Enzyme A and Enzyme E at room temperature (15 - 25°C) for immediate use or overnight at 2 - 8°C. Do not thaw at 37°C.

NOTE: Once thawed, use immediately or aliquot and store at -20°C until the expiry date as indicated on the label. After thawing the aliquots, use immediately. Do not re-freeze.

## SPLEEN DISSOCIATION ENZYME COCKTAIL

Prepare spleen dissociation enzyme cocktail fresh before use.

For one spleen, prepare 2.5 mL of spleen dissociation enzyme cocktail by combining the following in a STEMprep™ Sample Tube:

- 2.425 mL of Enzyme Diluent Z
- 50 µL of Enzyme A
- 25 µL of Enzyme E


## Sample Preparation

Dissect mouse spleen into a non-treated culture dish containing cold PBS. Remove any connective tissue before proceeding.

## Directions for Use

1. Transfer the spleen tissue (70 - 140 mg) into the STEMprep™ Sample Tube containing 2.5 mL of spleen dissociation enzyme cocktail.  
NOTE: Plastic Rack for Centrifuge Tubes, 50 mL (Catalog #200-0651) is provided with each STEMprep™ Tissue Dissociator System to manage STEMprep™ Sample Tubes while loading or removing tissues.
2. Tightly close the STEMprep™ Sample Tube lid (ensuring proper closure) and insert it into a slot on the STEMprep™ Tissue Dissociator. For complete instructions on using the instrument, refer to the Technical Manual: STEMprep™ Tissue Dissociator (Document #10000030598), available at [www.stemcell.com](http://www.stemcell.com), or contact us to request a copy.  
NOTE: For best results, avoid pinning the tissue with the rotor in the center of the STEMprep™ Sample Tube. Ensure total sample volume containing tissue does not exceed the "MAX" fill line on the tube.
3. Select and run the "Mouse Spleen" dissociation protocol on the STEMprep™ Tissue Dissociator.  
NOTE: If the tube will not be removed immediately after protocol completion, add a 4°C hold to preserve sample integrity.
4. After the protocol is complete, remove the STEMprep™ Sample Tube from the slot. Place a 70 µm cell strainer in a 50 mL conical tube. Pre-wet the strainer with 5 mL of EasySep™ Buffer. Carefully pour the sample over the strainer.
5. Rinse the STEMprep™ Sample Tube with 10 mL of cold EasySep™ Buffer and pour the wash over the strainer. Repeat this wash step once more. Top up the conical tube to 50 mL with EasySep™ Buffer and screw on the cap. Discard the strainer and STEMprep™ Sample Tube.  
NOTE: Small tissue fragments may remain in the filter after dissociation; these typically do not affect yield. Larger tissue pieces can be gently pushed through the strainer using the rubber end of a syringe plunger.
6. Centrifuge the conical tube at 300 x g for 10 minutes at room temperature with the brake on low. After centrifugation, carefully aspirate the supernatant.
7. OPTIONAL: If desired, perform red blood cell lysis by adding 5 mL of Ammonium Chloride Solution (Catalog #07850) to the cell pellet. Thoroughly mix the cell suspension by pipetting up and down. Incubate on ice for 10 minutes. Top up sample to 50 mL with EasySep™ Buffer. Centrifuge at 300 x g for 10 minutes at room temperature with the brake on low. After centrifugation, carefully aspirate the supernatant.  
NOTE: If proceeding to EasySep™, do not perform red blood cell lysis.
8. Resuspend cell pellet in the desired volume of EasySep™ Buffer or medium of your choice. Cells are ready for downstream use.

## Data

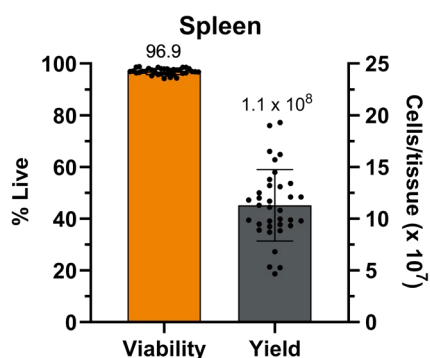

**Figure 1. Cell Viability and Yield Obtained Using the STEMprep™ Mouse Spleen Dissociation Kit**

Spleens from healthy C57BL/6 mice were dissociated into single-cell suspensions using the STEMprep™ Mouse Spleen Dissociation Kit on the STEMprep™ Tissue Dissociator. Following STEMprep™ processing, viability and yield of single-cell suspensions were assessed by flow cytometry. Data are presented as mean ± SD.

## Related Products

For more information about STEMprep™ kits and protocols, visit [www.stemcell.com/stemprep](http://www.stemcell.com/stemprep), contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED. FOR ADDITIONAL INFORMATION ON QUALITY AT STEMCELL, REFER TO [WWW.STEMCELL.COM/COMPLIANCE](http://WWW.STEMCELL.COM/COMPLIANCE).

Copyright © 2025 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, and STEMprep are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.