

Human Pluripotent Stem Cell Naïve State qPCR Array

For characterization of gene expression associated with hPSC naïve or primed states



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

The Human Pluripotent Stem Cell (hPSC) Naïve State Quantitative Polymerase Chain Reaction (qPCR) Array is designed for characterization of hPSCs and their status in the spectrum from naïve to primed pluripotency. Naïve state hPSCs are self-renewing, retain the characteristics of the pre-implantation blastocyst, have the capacity to differentiate to all somatic lineages without lineage bias, and demonstrate potential to differentiate into the germ lineage. Primed hPSCs are also self-renewing and have the ability to contribute to all somatic lineages, however some primed hPSC lines show lineage-specific differentiation bias. The hPSC Naïve State qPCR Array is designed for characterization of the gene expression profile of primed and naïve state hPSCs. Genes were selected based on their demonstrated differential expression in primed and naïve state hPSCs (Chan et al.; Davidson et al.; Gafni et al.; Theunissen et al.; Takashima et al.) or in hPSC-derived early ectodermal, endodermal, and mesodermal lineage cells (Huang et al.).

qPCR is a method for determining changes in steady-state mRNA levels of gene expression across multiple samples, generally normalized to the relative expression of internal control genes. Gene-specific primers are used in PCR to amplify target sequences within cDNA pools reverse-transcribed from mRNA. These PCR products contain hybridized sequence-specific probes that provide a fluorescent signal. Similar to TaqMan® technology, the fluorescent signal results from the 5' exonuclease activity of the Taq DNA polymerase on the probe, which is labeled with a reporter fluorophore at the 5' end and a quencher fluorophore at the 3' end. The rate of accumulation of the fluorescent signal is used to quantify the amount of cDNA present in the sample, and thereby the amount of mRNA present in the original cell lysate.

This qPCR Array contains validated primers and probes for detection of 90 genes whose expression is correlated with naïve and primed PSCs and their derivatives. There are also 6 wells containing primers and probes for endogenous (housekeeping) control genes. TATA box-binding protein (TBP) qPCR Array Control Template is provided separately as a synthetic DNA positive control, for use in a control well containing primers and probes for TBP.

An annotated list of genes, as well as plate layouts and software for analysis of qPCR results, are available at www.stemcell.com/qPCRanalysis.

Ordering Information

All kits listed below include TBP qPCR Control Template (Component #07518). For instrument compatibility, visit www.stemcell.com/naiveqPCRinstruments.

KIT CATALOG #	PLATE COMPONENT #	SIZE
07521	07503.1	1 Plate (96 wells)
07522	07503.2	1 Plate (96 wells)
07523	07503.3	1 Plate (96 wells)
07524	07503.4	1 Plate (96 wells)

Storage and Stability

Store plates at -20°C. Stable until expiry date (EXP) on boxtop label.

Store TBP qPCR Control Template at -20°C. Stable until expiry date (EXP) on label.

NOTE: Components may be shipped at room temperature (15 - 25°C) but should be stored at -20°C as indicated above.

Materials Required But Not Included

PRODUCT NAME	CATALOG #
qPCR Master Mix Kit - qPCR Master Mix (1 mL or 5 mL) - ROX Reference Dye (200 μL)	07516 (1 mL kit) OR 07517 (5 mL kit)
STEMscript™ cDNA Synthesis Kit with Oligo(dT) Primers OR STEMscript™ cDNA Synthesis Kit with Random Primers	79003 OR 79004
Nuclease-Free Water (not DEPC-treated)	79001
Optical adhesive film	e.g. Thermo Fisher 4311971

Directions for Use

Please read the entire protocol before proceeding.

Use sterile techniques when performing the following protocols.

Isolate RNA using standard protocols. Quantify RNA by optical density at 260 nm, determine purity using $A_{260/280}$, then convert to cDNA using a STEMscript™ cDNA Synthesis Kit.

Store cDNA at -20°C.

NOTE: Optimal concentration of cDNA for qPCR amplification is 20 - 100 ng/ μ L.

A. PREPARATION OF TBP qPCR CONTROL TEMPLATE AND cDNA COCKTAIL

- Thaw qPCR Master Mix, cDNA, and ROX Reference Dye (if using) on ice.
- If using ROX Reference Dye, add to qPCR Master Mix according to Table 1. For instruments not listed, refer to the manufacturer's instructions.

Table 1. Volume of ROX Reference Dye to add to qPCR Master Mix

PCR SYSTEM	VOLUME OF ROX REFERENCE DYE (μ L)	
	1 mL qPCR Master Mix	5 mL qPCR Master Mix
7900HT Fast (Applied Biosystems)	40	200
ViiA™ 7 (Applied Biosystems)	4	20
QuantStudio™ Flex (Applied Biosystems)	4	20

- Swirl bottle of qPCR Master Mix to mix thoroughly.
- Prepare **TBP qPCR Control Template** as follows:
 - Centrifuge TBP qPCR Control Template at 3000 x *g* for 3 - 5 seconds to pellet material to the bottom of the vial.
 - Add 20 μ L of nuclease-free water to the vial. Vortex the vial gently and thoroughly to resuspend the pellet.
 - Centrifuge at 3000 x *g* for 3 - 5 seconds to bring the liquid to the bottom of the vial.
- Prepare **cDNA Cocktail** as follows:
 - Mix cDNA by gently pipetting up and down. Centrifuge at 3000 x *g* for 3 - 5 seconds to bring liquid to the bottom of the vial.
 - To a 15 mL conical tube (e.g. Catalog #38009), add components according to Table 2.

Table 2. Preparation of cDNA Cocktail

cDNA Cocktail Components	VOLUME (μ L)	
	1 well	96 wells*
cDNA	1	108
qPCR Master Mix	5	540
Nuclease-free water	4	432
Total Volume	10	1080

*12.5% excess volume added to account for pipetting dead volume

- c. Cap the tube then gently vortex to mix thoroughly.
- d. Centrifuge at 3000 x *g* for 3 - 5 seconds to bring the liquid to the bottom of the tube.

B. PREPARATION OF qPCR PLATE

1. Carefully remove qPCR array plate from the box and plastic bag. Leave adhesive seal attached.
2. Centrifuge the plate at 1000 x *g* for 1 minute in a swinging bucket rotor fitted with plate holders.
3. Carefully remove and discard the adhesive seal on the plate.
4. Using a multichannel pipettor (e.g. Catalog #38064) and reagent reservoir (e.g. Catalog #38080), dispense reagents (from section A) into the plate wells as described below.
 - 5 μ L **TBP qPCR Control Template** + 5 μ L **qPCR Master Mix** in well H12 (see Figure 1)
 - 10 μ L **cDNA Cocktail** in all other wells

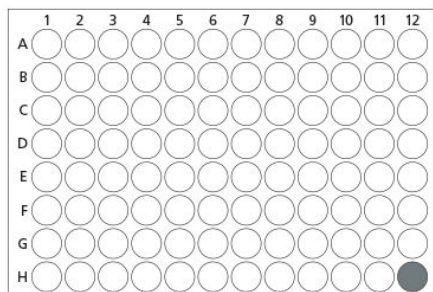


Figure 1. 96-Well Plate Diagram Indicating Well Containing TBP qPCR Control Template

5. Carefully cover and seal the plate using optical adhesive film.
6. Centrifuge the sealed plate at 1000 x *g* for 1 minute at room temperature (15 - 25°C) to remove bubbles from the bottom of the wells.
NOTE: Bubbles in the bottom of the wells will interfere with results.
7. Place the plate on ice.

C. qPCR

1. Program the thermocycler according to manufacturer's instructions.
2. If ROX Reference Dye is being used, calibrate thermocycler.
3. Add plate and run PCR program.
4. Save file including Ct (cycle threshold) values.
5. Import the Ct data from the qPCR instrument to the analysis tool available at www.stemcell.com/qPCRanalysis. This analysis tool can rapidly and accurately quantitate relative gene expression, and the user can change analysis settings with ease.

Related Products

For related products, including specialized cell culture and storage media, matrices, antibodies, cytokines, and small molecules, visit www.stemcell.com/hPSCworkflow or contact us at techsupport@stemcell.com.

References

- Chan Y-S et al. (2013) Induction of a human pluripotent state with distinct regulatory circuitry that resembles preimplantation epiblast. *Cell Stem Cell* 13(6): 663–75.
- Davidson KC et al. (2015) The pluripotent state in mouse and human. *Development* 142(18): 3090–9.
- Gafni O et al. (2013) Derivation of novel human ground state naive pluripotent stem cells. *Nature* 504(7479): 282–6.
- Huang G et al. (2015) Molecular basis of embryonic stem cell self-renewal: from signaling pathways to pluripotency network. *Cell Mol Life Sci* 72(9): 1741–57.
- Takashima Y et al. (2014) Resetting transcription factor control circuitry toward ground-state pluripotency in human. *Cell* 158(6): 1254–69.
- Theunissen TW et al. (2014) Systematic identification of culture conditions for induction and maintenance of naive human pluripotency. *Cell Stem Cell* 15(4): 471–87.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2018 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. TaqMan is a registered trademark of Roche Molecular Systems, Inc. QuantStudio and ViiA are trademarks of Thermo Fisher Scientific Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.