

Anti-Mouse EPCR Antibody, Clone RMEPCR1560 (1560), FITC

Rat monoclonal IgG2b antibody against mouse EPCR (CD201), FITC-conjugated

Catalog #60038Fl.1	25 Tests	20 µL/test
Catalog #100-1576	100 Tests	20 µL/test

Product Description

This monoclonal antibody reacts with the endothelial protein C receptor (EPCR or CD201), an ~25 kDa type I transmembrane glycoprotein expressed by endothelial cells, subsets of hematopoietic stem cells (HSCs) and dendritic cells, and several malignant cell lines. It is also found in a soluble form in plasma. EPCR exhibits homology with the MHC class 1/CD1 protein family. EPCR binds protein C and activated protein C, thus augmenting protein C activation by the thrombin-thrombomodulin complex and regulating blood coagulation and inflammation. EPCR protein expression has been detected on ~1.5% of mouse bone marrow cells. Purified EPCR+ cells are highly enriched for HSC activity, as evidenced by high in vivo repopulation activity. Moreover, EPCR expression is associated with the stem cell activity of bone marrow cell populations isolated using conventional markers, indicating the usefulness of EPCR as a single marker for the identification of mouse HSCs.

Target Antigen:	EPCR (CD201)
Alternative Names:	Activated protein C receptor, APC receptor, CCD41, CD201, endothelial protein C receptor, PROCR, protein C receptor
Gene ID:	19124
Species Reactivity:	Mouse
Host Species:	Rat
Clonality:	Monoclonal
Clone:	RMEPCR1560 (1560)
Isotype:	IgG2b, kappa
Immunogen:	Soluble form of mouse EPCR protein
Conjugate:	FITC (Fluorescein isothiocyanate)

Applications

Verified Applications: FC

Reported Applications: FC

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FACS: Fluorescence-activated cell sorting; FC: Flow cytometry; FCXM: Flow cytometric crossmatch assay; FISH: Fluorescence in situ hybridization; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IHC-F: Immunohistochemistry (frozen-tissue); IHC-P: Immunohistochemistry (paraffin-embedded); IP: Immunoprecipitation; NMR: Nuclear magnetic resonance spectroscopy; RIA: Radioimmunoassay; WB: Western blotting

Properties


Product Formulation: Phosphate-buffered saline containing less than 0.1% (w/v) sodium azide and less than 0.1% (w/v) bovine serum albumin

Purification: The antibody was purified by affinity chromatography and conjugated with FITC under optimal conditions.

Stability and Storage: Product stable at 2 - 8°C when stored undiluted. Do not freeze. Protect product from prolonged exposure to light. For product expiry date, contact techsupport@stemcell.com.

Directions for Use: For flow cytometry, the suggested use of this antibody is 20 µL per 1×10^6 cells in 100 µL or per 100 µL of whole blood. It is recommended that the antibody be titrated for optimal performance for each application. Not recommended for use with mouse bone marrow cells.

Data

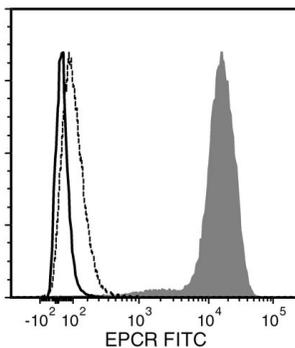

Flow cytometry analysis of HEK-293 mEPCR-transfected cells (filled histogram) or non-transfected HEK-293 cells (negative control cells, dashed line histogram), labeled with Anti-Mouse EPCR Antibody, Clone RMEPCR1560 (1560), FITC. Labeling of HEK-293 mEPCR-transfected cells with a rat IgG2b, kappa FITC isotype control antibody is shown (solid line histogram).

Related Products

For a complete list of antibodies, including other conjugates, sizes, and clones, as well as related products available from STEMCELL Technologies, visit www.stemcell.com/antibodies, or contact us at techsupport@stemcell.com.

References

- Akashi K et al. (2003) Transcriptional accessibility for genes of multiple tissues and hematopoietic lineages is hierarchically controlled during early hematopoiesis. *Blood* 101(2): 383–9.
- Balazs AB et al. (2006) Endothelial protein C receptor (CD201) explicitly identifies hematopoietic stem cells in murine bone marrow. *Blood* 107(6): 2317–21.
- Crawley JTB et al. (2002) Distribution of endothelial cell protein C/activated protein C receptor (EPCR) during mouse embryo development. *Thromb Haemost* 88(2): 259–66.
- Esmon CT. (2001) The normal role of activated protein C in maintaining homeostasis and its relevance to critical illness. *Crit Care* 5(2): S7-12.
- Gu J-M et al. (2002) Disruption of the endothelial cell protein C receptor gene in mice causes placental thrombosis and early embryonic lethality. *J Biol Chem* 277(45): 43335–43.
- Kent DG et al. (2013) Self-renewal of single mouse hematopoietic stem cells is reduced by JAK2V617F without compromising progenitor cell expansion. *PLoS Biol* 11(6): e1001576.
- Kirschner K et al. (2017) Proliferation drives aging-related functional decline in a subpopulation of the hematopoietic stem cell compartment. *Cell Rep* 19(8): 1503–11.
- Liaw PC et al. (2000) Mechanisms by which soluble endothelial cell protein C receptor modulates protein C and activated protein C function. *J Biol Chem* 275(8): 5447–52.
- Mayle A et al. (2013) Flow cytometry analysis of murine hematopoietic stem cells. *Cytom Part A* 83A(1): 27–37. (FACS, FC) 5. Gasparetto M et al. (2012) Aldehyde dehydrogenases are regulators of hematopoietic stem cell numbers and B-cell development. *Exp Hematol* 40(4): 318–29.e2.
- Ramalho-Santos M et al. (2002) “Stemness”: transcriptional profiling of embryonic and adult stem cells. *Science* 298(5593): 597–600.
- Sundaram J et al. (2014) Factor VIIa binding to endothelial cell protein C receptor protects vascular barrier integrity in vivo. *J Thromb Haemost* 12 (5): 690–700.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2024 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.