

Mouse Recombinant M-CSF (E.coli-expressed)

Macrophage colony-stimulating factor

Catalog #78059.1	10 µg
Catalog #78059	100 µg
Catalog #78059.3	500 µg
Catalog #78059.2	1000 µg

Product Description

Macrophage colony-stimulating factor (M-CSF) is a homodimeric glycoprotein growth factor that regulates proliferation and differentiation of myeloid hematopoietic progenitors to mononuclear phagocytic cell lineages, including monocytes, macrophages, and osteoclasts. M-CSF is a crucial factor for the development of tissue-resident macrophages in most tissues (Ginhoux & Jung). It is required for the maturation and activation of monocytes and macrophages, and regulates inflammatory responses in conjunction with other stimuli such as IFN- γ , LPS, and IL-4 (Murray et al.). M-CSF is also required for bone resorption by osteoclasts, and is involved in the development and regulation of placenta, mammary gland, and brain. M-CSF is produced by monocytes, fibroblasts, osteoclasts, stromal cells, endothelial cells, and tumor cells (Chockalingam & Ghosh).

M-CSF exerts its biological effects by signaling through a receptor tyrosine kinase (CSF-1R or M-CSF-R) encoded by the *c-fms* proto-oncogene (Hamilton). CSF-1R shares similar structural features with other growth factor receptors, including the stem cell factor (SCF) receptor, platelet-derived growth factor receptor (PDGF-R), and Flt3/Flk-2 receptor tyrosine kinase. Stimulation of the CSF-1R upon binding to M-CSF activates MAPK, PI3K, and PLC γ signaling pathways (Chockalingam & Ghosh).

Product Information

Alternative Names:	Colony stimulating factor 1, CSF-1
Accession Number:	P07141
Amino Acid Sequence:	MKEVSEHCSH MIGNGHLKVL QQLIDSQMET SCQIAFEFVD QEQLDDPVCY LKKAFFLVQD IIDETMRFKD NTPNANATER LQELSNNLNS CFTKDYEEQN KACVRTFHET PLQLEKIKN FFNETKNLLE KDOWNIFTKNC NNSFAKCSSR DVVTKP
Predicted Molecular Mass:	18.2 kDa monomer; 36.4 kDa dimer
Species:	Mouse
Product Formulation:	Lyophilized from a sterile-filtered aqueous solution containing sodium phosphate and sodium chloride, pH 7.5.
Source:	E. coli
Purity:	≥ 95%


Specifications

Activity:	The specific activity is $\geq 1 \times 10^5$ units/mg ($EC_{50} \leq 10$ ng/mL), as determined by a cell proliferation assay using NFS-60 cells.
Endotoxin Level:	Measured by kinetic Limulus amebocyte lysate (LAL) analysis and is ≤ 1 EU/ μ g protein.

Preparation and Storage

Stability and Storage:	Store at -20 to -80°C. Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Reconstitute the product in sterile water to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. OPTIONAL: After reconstitution, if product will not be used immediately, dilute with concentrated bovine serum albumin (BSA) to a final BSA concentration of 0.1%. The effect of storage of stock solution on product performance should be tested for each application. As a general guide, do not store at 2 - 8°C for more than 1 month or at -20 to -80°C for more than 3 months. Avoid repeated freeze-thaw cycles.

Data

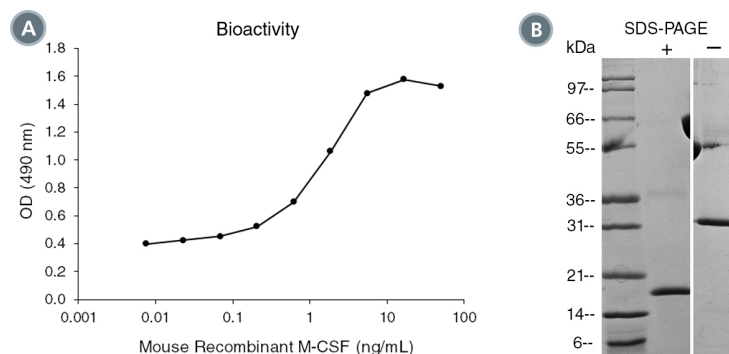

(A) The biological activity of Mouse Recombinant M-CSF was tested by its ability to promote the proliferation of NFS-60 cells. Cell proliferation was measured after 44 hours of culture. The EC50 is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC50 in the above example is 1.1 - 1.6 ng/mL. (B) 1 µg of Mouse Recombinant M-CSF was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Mouse Recombinant M-CSF is a homodimer of 18.2 kDa subunits with a predicted total predicted molecular mass of 36.4 kDa.

Related Products

For a complete list of cytokines or peptide pools, as well as related products available from STEMCELL Technologies, visit www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Chockalingam S & Ghosh SS. (2014) Macrophage colony-stimulating factor and cancer: a review. *Tumour Biol* 35(11): 10635-44.
- Ginhoux F & Jung S. (2014) Monocytes and macrophages: developmental pathways and tissue homeostasis. *Nat Rev Immunol* 14(6): 392-404.
- Hamilton JA. (1997) CSF-1 signal transduction. *J Leukoc Biol* 62(2): 145-55.
- Murray PJ et al. (2014) Macrophage Activation and Polarization: Nomenclature and Experimental Guidelines. *Immunity* 41(1): 14-20.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2024 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.