

# Human Recombinant IL-3, ACF

Interleukin 3, animal component-free

Catalog #78146	10 µg
Catalog #78146.1	100 µg
Catalog #78146.2	1000 µg

## Product Description

Interleukin 3 (IL-3) is a species-specific pleiotropic cytokine that promotes the survival and proliferation of pluripotent hematopoietic stem cells and lineage-committed progenitor cells and their differentiation into mature cells of most lineages, including basophils, neutrophils, eosinophils, macrophages, dendritic cells, erythrocytes, and megakaryocytes (Broughton et al.; Dorssers et al.; Yang et al.). IL-3 is produced by activated T cells and has a physiological role in inflammation and allergies by promoting the secretion of inflammatory mediators such as histamine, IL-4, and IL-6 by basophils and eosinophils (Broughton et al.). The IL-3 receptor consists of a unique alpha subunit (CD123) and a beta common subunit (βc or CD131) that is shared with the receptors for IL-5 and GM-CSF, and is the principal signal transduction subunit for these cytokines. IL-3 binding to the heterodimeric receptor activates JAK/STAT, MAPK, and PI3K signaling pathways (Woodcock et al.). This product is animal component-free.

## Product Information

Alternative Names:	FGFA, Fibroblast growth factor-10, Keratinocyte growth factor-2
Accession Number:	P08700
Amino Acid Sequence:	MAPMTQTTP L KTSWVNCSNM IDEIITHLKQ PPLPLDFNN LNGEDQDILM ENNLRRPNLE AFNRAVKSLQ NASAIESILK NLLPCLPLAT AAPTRHPIHI KDGDWNEFRR KLTFYKLTLE NAQAQQTTL S LAIF
Predicted Molecular Mass:	15.2 kDa
Species:	Human
Product Formulation:	Lyophilized from a sterile-filtered aqueous solution containing sodium phosphate and sodium chloride, pH 7.5.
Source:	E. coli
Purity:	≥ 97%


## Specifications

Activity:	The EC <sub>50</sub> is ≤ 2 ng/mL, as determined by a cell proliferation assay using TF-1 cells. The specific activity is approximately 4.2 x 10 <sup>3</sup> IU/μg, as calibrated against the human recombinant IL-3 WHO International Standard (NIBSC code: 91/510).
Endotoxin Level:	Measured by kinetic Limulus amoebocyte lysate (LAL) analysis and is ≤ 0.1 EU/μg protein.

## Preparation and Storage

Stability and Storage:	Store at -20 to -80°C. Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Reconstitute the product in sterile water to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. OPTIONAL: After reconstitution, if product will not be used immediately, dilute with concentrated bovine serum albumin (BSA) to a final BSA concentration of 0.1%. The effect of storage of stock solution on product performance should be tested for each application. As a general guide, do not store at 2 - 8°C for more than 1 month or at -20 to -80°C for more than 3 months. Avoid repeated freeze-thaw cycles.

## Data

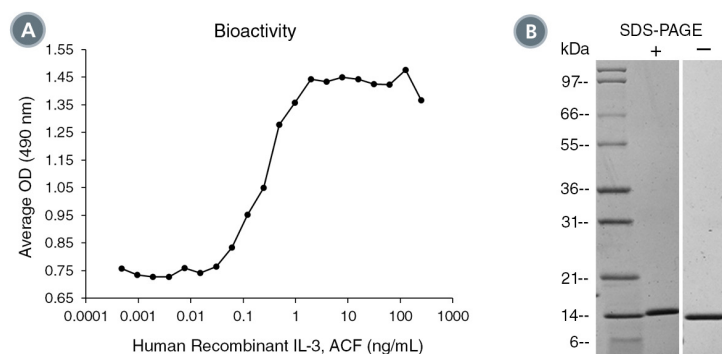

**Figure 1. Biological Activity and Molecular Mass of Human Recombinant IL-3, ACF**

(A) The biological activity of Human Recombinant IL-3, ACF was tested by its ability to promote the proliferation of TF-1 cells. The EC<sub>50</sub> is defined as the effective concentration of the cytokine at which cell proliferation is at 50% of maximum. The EC<sub>50</sub> in the example above is 0.244 ng/mL. (B) 1 µg of Human Recombinant IL-3, ACF was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Human Recombinant IL-3, ACF has a predicted molecular mass of 15.2 kDa.

## Related Products

For a complete list of cytokines or peptide pools, as well as related products available from STEMCELL Technologies, visit [www.stemcell.com/cytokines](http://www.stemcell.com/cytokines) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

## References

- Broughton SE et al. (2012) The GM-CSF/IL-3/IL-5 cytokine receptor family: from ligand recognition to initiation of signaling. *Immunol Rev* 250(1): 277–302.
- Dorssers L et al. (1987) Characterization of a human multilineage-colony-stimulating factor cDNA clone identified by a conserved noncoding sequence in mouse interleukin-3. *Gene* 55(1): 115–24.
- Woodcock JM et al. (1999) The functional basis of granulocyte-macrophage colony stimulating factor, interleukin-3 and interleukin-5 receptor activation, basic and clinical implications. *Int J Biochem Cell Biol* 31(10): 1017–25.
- Yang YC et al. (1986) Human IL-3 (multi-CSF): identification by expression cloning of a novel hematopoietic growth factor related to murine IL-3. *Cell* 47(1): 3–10.

PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2024 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.