

**Anti-Human CD105 Antibody,
Clone 43A3, Alexa Fluor® 488**



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Antibodies

Mouse monoclonal IgG1 antibody
against human CD105 (endoglin),
Alexa Fluor® 488-conjugated

Catalog #60039AD
#60039AD.1

100 tests	5 µL/test
25 tests	5 µL/test

FOR RESEARCH USE ONLY. NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES.

Product Description

The 43A3 antibody reacts with CD105 (endoglin), an ~180 kDa cell surface glycoprotein which is a disulfide-bonded homodimer of ~90 kDa type I transmembrane subunits. CD105 is a component of the TGF- β receptor complex and is expressed by vascular endothelial smooth muscle cells, syncytiotrophoblasts of placenta and activated macrophages, and at relatively low levels by stromal fibroblasts. Its expression is also observed in some types of tumors and levels are up-regulated on the endothelium during angiogenesis. In concert with signaling receptors, CD105 binds to TGF- β 1 and TGF- β 3 with high affinity, but does not bind TGF- β 2. Other ligands reportedly include Activin A, BMP-2 and BMP-7. CD105 has important roles in angiogenesis, cardiovascular development and vascular remodeling, and the protein serves a regulatory role in cytoskeletal reorganization by modulating the sites of focal adhesion and cellular migration. Certain mutations in CD105 result in the autosomal dominant disorder hereditary hemorrhagic telangiectasia.

Target Antigen Name:	CD105 (Endoglin)
Alternative Names:	Endoglin
Gene ID:	2022
Species Reactivity:	Human, Mouse
Host Species:	Mouse
Clonality:	Monoclonal
Clone:	43A3
Isotype:	IgG1, kappa
Immunogen:	L-cells transfected with human CD105
Conjugate:	Alexa Fluor® 488

Applications


Verified:	FC
Reported:	FC
Special Applications:	This antibody clone has been verified for labeling human mesenchymal cells grown in MesenCult™ Proliferation Kit (Human; Catalog #05411) and MesenCult™-XF Medium (Catalog #05420).

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FC: Flow cytometry; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IP: Immunoprecipitation; WB: Western blotting

Properties

Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) bovine serum albumin (origin USA)
Purification:	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor® 488 under optimal conditions. The solution is free of unconjugated Alexa Fluor® 488.
Stability and Storage:	Product stable at 2 - 8°C when stored undiluted. Do not freeze. Protect product from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com .
Directions for Use:	For flow cytometry the suggested use of this antibody is 5 µL per 1 x 10 ⁶ cells in 100 µL volume or per 100 µL of whole blood. It is recommended that the antibody be titrated for optimal performance for each application.

Data

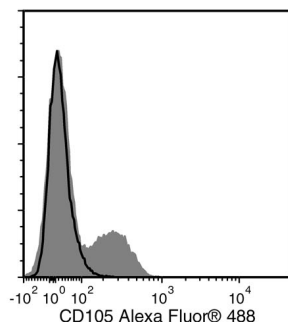

Flow cytometry analysis of human peripheral blood mononuclear cells (PBMCs) labeled with Anti-Human CD105 Antibody, Clone 43A3, Alexa Fluor® 488 (filled histogram) or a mouse IgG1, kappa Alexa Fluor® 488 isotype control antibody (solid line histogram).

Related Products

For a complete list of antibodies, including other conjugates, sizes and clones, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/antibodies or contact us at techsupport@stemcell.com.

References

1. Bühring HJ, et al. Endoglin is expressed on a subpopulation of immature erythroid cells of normal human bone marrow. *Leukemia* 5(10): 841-47, 1991
2. Wu X, et al. Cloning and characterization of the murine activin receptor like kinase-1 (ALK-1) homolog. *Biochem Biophys Res Commun* 216(1): 78-83, 1995
3. Altomonte M, et al. Expression and structural features of Endoglin (CD105), a transforming growth factor β 1 and β 3 binding protein, in human melanoma. *Br J Cancer* 74(10): 1586-91, 1996
4. Gallione CJ, et al. Mutation and expression analysis of the Endoglin gene in hereditary hemorrhagic telangiectasia reveals null alleles. *Hum Mutat* 11(4): 286-94, 1998
5. Azuma, H. Genetic and molecular pathogenesis of hereditary hemorrhagic telangiectasia. *J Med Invest* 47(3-4): 81-90, 2000
6. Bourdeau, A, et al. Endoglin-deficient mice, a unique model to study hereditary hemorrhagic telangiectasia. *Trends Cardiovasc Med* 10(7): 279-85, 2000
7. Pierelli L, et al. CD105 (endoglin) expression on hematopoietic stem/progenitor cells. *Leuk Lymphoma* 42(6): 1195-206, 2001
8. Mason D, et al. Eds. *Leucocyte Typing VII*. Oxford University Press. New York, 2002 (FC)
9. Vogel W, et al. Heterogeneity among human bone marrow-derived mesenchymal stem cells and neural progenitor cells. *Haematologica* 88(2): 126-33, 2002 (FC, ICC)
10. HonMou O, et al. Intravenous administration of auto serum-expanded autologous mesenchymal stem cells in stroke. *Brain* 134(6): 1790-807, 2011
11. Osaka M, et al. Intravenous administration of mesenchymal stem cells derived from bone marrow after contusive spinal cord injury improves functional outcome. *Brain Res* 1343: 226-35, 2010
12. Akiyama K, et al. Characterization of bone marrow derived mesenchymal stem cells in suspension. *Stem Cell Res Ther* 3: 40, 2012 (WB)
13. Marolt D, et al. Engineering bone tissue from human embryonic stem cells. *Proc Natl Acad Sci USA* 109(22): 8705-09, 2012 (IHC)

Copyright © 2013 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists and MesenCult are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. Alexa Fluor® is a registered trademark of Life Technologies Corporation. This product is licensed for internal research use only and its sale is expressly conditioned on the buyer not using it for manufacturing, performing a service, or medical test, or otherwise generating revenue. For use other than research, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485 MEDICAL DEVICE STANDARDS.