

Anti-Human CD34 Antibody, Clone 581

Antibodies

Mouse monoclonal IgG1 antibody
against human CD34, unconjugated



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #60013

100 µg 0.5 mg/mL

FOR RESEARCH USE ONLY. NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES.

Product Description

The 581 antibody reacts with human CD34, an ~105 - 120 kDa type 1 transmembrane glycoprotein expressed on the surface of most human hematopoietic stem and progenitor cells (HSPCs) as well as on mesenchymal stem cells, embryonic fibroblasts, endothelial cells, neurons, and some tumor cell lines. CD34 is expressed only transiently during hematopoiesis, so the frequency of CD34+ cells is low in bone marrow or cord blood (~1 - 5%) and very low (~0.1 - 0.5%) in peripheral blood. CD34 is a marker used to identify and isolate HSPCs capable of cell engraftment. CD34 is thought to mediate attachment of stem cells to the bone marrow extracellular matrix or directly to stromal cells during early hematopoiesis, and to be involved in lymphocyte recruitment through binding to the ligands L- and E-selectin. Distinct epitope groups have been assigned to CD34 based on their sensitivity to enzymatic cleavage, with the 581 antibody recognizing a class III epitope (resistant to neuraminidase and O-glycoprotease).

Target Antigen Name:	CD34
Alternative Names:	Gp105-120, My10
Gene ID:	947
Species Reactivity:	Human
Host Species:	Mouse
Clonality:	Monoclonal
Clone:	581
Isotype:	IgG1, kappa
Immunogen:	Human CD34+ leukemic cells
Conjugate:	Unconjugated

Applications


Verified:	FC
Reported:	CyTOF®, FC, ICC, IF, IHC
Special Applications:	This antibody clone has been verified for purity assessments of cells isolated with EasySep™ Human CD34 Positive Selection Kit (Catalog #18056) and for labeling human mesenchymal cells grown in MesenCult™ Proliferation Kit (Human; Catalog #05411).

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FC: Flow cytometry; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IP: Immunoprecipitation; WB: Western blotting

Properties

Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Purification:	The antibody was purified by affinity chromatography.
Stability and Storage:	Product stable at 2 - 8°C when stored undiluted. Do not freeze. For product expiry date, please request a lot-specific Certificate of Analysis from techsupport@stemcell.com .
Directions for Use:	For flow cytometry the suggested use of this antibody is ≤ 2.0 µg per 1 x 10 ⁶ cells in 100 µL volume. It is recommended that the antibody be titrated for optimal performance for each application. Clone 581 is not recommended for WB.

Data

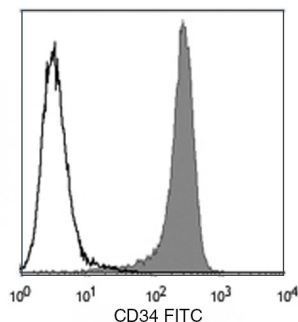

Flow cytometry analysis of human KG1a myeloid cells labeled with Anti-Human CD34 Antibody, Clone 581, followed by anti-mouse IgG, FITC (filled histogram) or a mouse IgG1, kappa isotype control antibody followed by anti-mouse IgG, FITC (open histogram).

Related Products

PRODUCT NAME	CATALOG #	SIZE
Anti-Human CD34 Antibody, Clone 581	60013	100 µg
Anti-Human CD34 Antibody, Clone 581, PE	60013PE	100 tests
Anti-Human CD34 Antibody, Clone 581, PE	60013PE.1	25 tests
Anti-Human CD34 Antibody, Clone 581, Alexa Fluor® 488	60013AD	100 tests
Anti-Human CD34 Antibody, Clone 581, Alexa Fluor® 488	60013AD.1	25 tests

References

1. Civin CI, et al. Antigenic analysis of hematopoiesis. III. A hematopoietic progenitor cell surface antigen defined by a monoclonal antibody raised against KG-1a cells. *J Immunol* 133(1): 157-65, 1984
2. Strauss LC, et al. Antigenic analysis of hematopoiesis. V. Characterization of My-10 antigen expression by normal lymphohematopoietic progenitor cells. *Exp Hematol* 14(9): 878-86, 1986
3. Terstappen LW, et al. Sequential generations of hematopoietic colonies derived from single nonlineage-committed CD34+CD38- progenitor cells. *Blood* 77(6): 1218-27, 1991 (FC)
4. Andrews RG, et al. CD34+ marrow cells, devoid of T and B lymphocytes, reconstitute stable lymphopoiesis and myelopoiesis in lethally irradiated allogeneic baboons. *Blood* 80(7): 1693-701, 1992 (FC, IF)
5. Schlossman SF, et al. Eds. *Leukocyte Typing V: White Cell Differentiation Antigens*. Oxford University Press, New York, pp. 840-46, 1995
6. Rudin CE, et al. Concordant expression of class II and class III CD34 epitopes on haemopoietic cells in leukapheresis and cord blood samples with CD45/CD34 dual staining. *Br J Haematol* 97(2): 488-91, 1997 (FC)
7. Steen R, Egeland T. CD34 molecule epitope distribution on cells of haematopoietic origin. *Leuk Lymphoma* 30(1-2): 23-30, 1998
8. Yoshino N, et al. Upgrading of flow cytometric analysis for absolute counts, cytokines and other antigenic molecules of Cynomolgus monkeys (*Macaca fascicularis*) by using anti-human cross-reactive antibodies. *Exp Anim* 49(2): 97-110, 2000 (FC)
9. Joseph A. Expression of CD34 and L-selectin on human corneal keratocytes. *Invest Ophthalmol Vis Sci* 44(11): 4689-92, 2003 (IHC)
10. Maumus M, et al. Native human adipose stromal cells: localization, morphology and phenotype. *Int J Obes* 35(9): 1141-53, 2011 (FC, ICC, IHC)

Copyright © 2013 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists and EasySep are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. Alexa Fluor® is a registered trademark of Life Technologies Corporation. This product is licensed for internal research use only and its sale is expressly conditioned on the buyer not using it for manufacturing, performing a service, or medical test, or otherwise generating revenue. For use other than research, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485 MEDICAL DEVICE STANDARDS.