

Anti-Mouse CD49f Antibody, Clone GoH3



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Antibodies

Rat monoclonal IgG2a antibody
against human, mouse, rhesus CD49f
(integrin $\alpha 6$), unconjugated

Catalog #60037

100 μ g 0.5 mg/mL

Product Description

The GoH3 antibody reacts with CD49f (integrin $\alpha 6$), an ~150 kDa transmembrane glycoprotein that associates non-covalently with CD29 (integrin $\beta 1$) or CD104 (integrin $\beta 4$) to form the heterodimeric receptors VLA-6 and $\alpha 6\beta 4$, which bind the extracellular matrix protein laminin. CD49f is a disulfide-linked dimer comprising an ~120 kDa heavy chain and an ~30 kDa membrane-bound light chain. Splice variants exist, which affect the cytoplasmic domain of the protein. CD49f is expressed on the surface of T cells, monocytes, platelets, placental trophoblasts, epithelial and endothelial cells. It is involved in cell adhesion and regulating signaling pathways involved in a variety of processes, including the activation and proliferation of T cells, and the differentiation and maintenance of pluripotency of stem cells. CD49f is considered the most important marker for selecting mouse mammary stem and progenitor cells. The GoH3 antibody reacts with an extracellular epitope on CD49f and reportedly blocks integrin $\alpha 6$ function in vivo and binding of integrin $\alpha 6$ to laminin in vitro.

Target Antigen Name:	CD49f (Integrin $\alpha 6$)
Alternative Names:	$\alpha 6$ integrin, integrin $\alpha 6$, VLA-6 α chain
Gene ID:	16403/3655
Species Reactivity:	Human, Mouse, Rhesus, Cynomolgus, Baboon, Chimpanzee, Capuchin Monkey, Cat, Cow, Dog, Horse, Pig, Rabbit, Sheep
Host Species:	Rat (SD)
Clonality:	Monoclonal
Clone:	GoH3
Isotype:	IgG2a, kappa
Immunogen:	Mouse (BALB/c) mammary tumor cells
Conjugate:	Unconjugated

Applications


Verified:	FC
Reported:	FA, FC, ICC, IF, IHC, IP, RIA

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FC: Flow cytometry; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IP: Immunoprecipitation; RIA: Radioimmunoassay; WB: Western blotting

Properties

Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide
Purification:	The antibody was purified by affinity chromatography.
Stability and Storage:	Product stable at 2 - 8°C when stored undiluted. Do not freeze. For product expiry date, please contact techsupport@stemcell.com .
Directions for Use:	For flow cytometry the suggested use of this antibody is $\leq 2 \mu$ g per 1×10^6 cells in 100 μ L volume. It is recommended that the antibody be titrated for optimal performance for each application.

Data

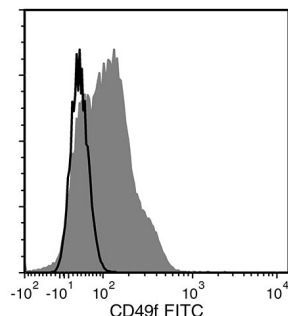

Flow cytometry analysis of human peripheral blood mononuclear cells (PBMCs) labeled with Anti-Mouse CD49f Antibody, Clone GoH3, followed by a mouse anti-rat IgG2a antibody, FITC (filled histogram), or Rat IgG2a, kappa Isotype Control Antibody, Clone RTK2758 (Catalog #60076), followed by a mouse anti-rat IgG2a antibody, FITC (solid line histogram).

Related Products

For a complete list of antibodies, including other conjugates, sizes and clones, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/antibodies or contact us at techsupport@stemcell.com.

References

1. Yu K-R et al. (2012) CD49f enhances multipotency and maintains stemness through the direct regulation of OCT4 and SOX2. *Stem Cells* 30(5): 876–87.
2. Grange C et al. (2011) Microvesicles released from human renal cancer stem cells stimulate angiogenesis and formation of lung premetastatic niche. *Cancer Res* 71(15): 5346–56. (FC)
3. Jo M et al. (2010) Cell signaling by urokinase-type plasminogen activator receptor induces stem cell-like properties in breast cancer cells. *Cancer Res* 70(21): 8948–58. (FC)
4. Stingl J et al. (2006) Purification and unique properties of mammary epithelial stem cells. *Nature* 439(7079): 993–7. (FC)
5. Ambrose HE & Wagner SD. (2004) Alpha6-integrin is expressed on germinal centre B cells and modifies growth of a B-cell line. *Immunology* 111(4): 400–6. (FC)
6. Yoshino N et al. (2000) Upgrading of flow cytometric analysis for absolute counts, cytokines and other antigenic molecules of cynomolgus monkeys (*Macaca fascicularis*) by using anti-human cross-reactive antibodies. *Exp Anim* 49(2): 97–110. (FC)
7. Lee EC et al. (1992) The integrin alpha 6 beta 4 is a laminin receptor. *J Cell Biol* 117(3): 671–8. (Blocking, FA, IP)
8. Aumailley M et al. (1990) Antibody to integrin alpha 6 subunit specifically inhibits cell-binding to laminin fragment 8. *Exp Cell Res* 188(1): 55–60. (FA)
9. Hemler ME et al. (1988) Multiple very late antigen (VLA) heterodimers on platelets. Evidence for distinct VLA-2, VLA-5 (fibronectin receptor), and VLA-6 structures. *J Biol Chem* 263(16): 7660–5.
10. Sonnenberg A et al. (1986) Development of mouse mammary gland: identification of stages in differentiation of luminal and myoepithelial cells using monoclonal antibodies and polyvalent antiserum against keratin. *J Histochem Cytochem* 34(8): 1037–46. (FC, IHC, IF)

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists and StemSpan are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. Alexa Fluor® is a registered trademark of Life Technologies Corporation. This product is licensed for internal research use only and its sale is expressly conditioned on the buyer not using it for manufacturing, performing a service, or medical test, or otherwise generating revenue. For use other than research, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.