



EasySep™ Mouse/Human Chimera Isolation Kit

Negative Selection

Catalog #19849

For processing 1×10^9 cells



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Document #28055 | Version 1_1_0

Description

Isolate untouched and highly purified human cells from mouse splenocytes, bone marrow, or peripheral blood by immunomagnetic negative selection. When using single-cell suspensions from other tissue types, this kit may require optimization.

- Fast, easy-to-use and column-free
- Isolated cells are untouched

This kit targets non-hematopoietic mouse cells for removal with biotinylated antibodies recognizing specific cell surface markers. Unwanted cells are labeled with biotinylated antibodies and streptavidin-coated magnetic particles, and separated without columns using an EasySep™ magnet. Desired cells are simply poured off into a new tube. Isolated cells are immediately available for downstream applications such as flow cytometry, culture or DNA/RNA extraction.

Component Descriptions

COMPONENT NAME	COMPONENT #	QUANTITY	STORAGE	SHELF LIFE	FORMAT
EasySep™ Mouse/Human Chimera Isolation Cocktail	19849C	1 x 0.5 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A combination of monoclonal antibodies in PBS and 0.1% BSA.
EasySep™ Streptavidin Rapidspheres™ 50001	50001	1 x 1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A suspension of magnetic particles in PBS.
Normal Rat Serum	13551	1 x 2 mL	Store at -20°C.	Stable until expiry date (EXP) on label.	Mycoplasma-free normal rat serum.

BSA - bovine serum albumin; PBS - phosphate-buffered saline

Components may be shipped at room temperature (15 - 25°C) but should be stored as indicated above.

Additional Reagent Stability Information

REAGENT NAME	STORAGE	SHELF LIFE
Normal Rat Serum (in-use)	Store at 2 - 8°C.	Stable for at least 2 months. Do not exceed expiry date (EXP) on label.

Sample Preparation

BONE MARROW

Flush bone marrow cells from femur and tibia into recommended medium using a syringe equipped with a 23 gauge needle. Disperse clumps by gently passing the cell suspension through the syringe several times. Alternatively, crush bones using a mortar and pestle. Remove remaining clumps and debris by passing cell suspension through a 70 µm mesh nylon strainer. Centrifuge at 300 x g for 10 minutes and resuspend cells at 1×10^8 cells/mL in recommended medium.

SPLEEN

Disrupt spleen in PBS or Hanks' Balanced Salt Solution (HBSS) containing 2% fetal bovine serum (FBS). Remove clumps and debris by passing cell suspension through a 70 µm mesh nylon strainer. Centrifuge at 300 x g for 10 minutes and resuspend at 1×10^8 nucleated cells/mL in recommended medium.

Ammonium chloride treatment is not recommended when preparing the cells for separation.

PERIPHERAL BLOOD

Blood should be lysed prior to use. Mix 1 part blood with 9 parts Ammonium Chloride Solution (Catalog #07800) and incubate on ice for 15 minutes. Centrifuge at 300 x g for 6 minutes. Discard supernatant and wash cell pellet once with recommended medium. Discard supernatant and resuspend cell pellet at 1×10^8 cells/mL in recommended medium.



Recommended Medium

EasySep™ Buffer (Catalog #20144), RoboSep™ Buffer (Catalog #20104), or PBS (Catalog #37350) containing 2% FBS and 1 mM EDTA. HBSS (Catalog #37250) can be used in place of PBS. Medium should be free of Ca^{++} , Mg^{++} , and biotin.

Directions for Use – Manual EasySep™ Protocols

See page 1 for Sample Preparation and Recommended Medium. Refer to Tables 1 and 2 for detailed instructions regarding the EasySep™ procedure for each magnet.




Table 1. EasySep™ Mouse/Human Chimera Isolation Kit Protocol

		EASYSEP™ MAGNETS	
STEP	INSTRUCTIONS	 EasySep™ (Catalog #18000)	 “The Big Easy” (Catalog #18001)
1	Prepare sample at the indicated cell concentration within the volume range.	1 x 10 ⁸ cells/mL 0.25 - 2 mL	1 x 10 ⁸ cells/mL 0.5 - 8 mL
2	Add Rat Serum to sample.	50 µL/mL of sample	50 µL/mL of sample
3	Add sample to required tube.	5 mL (12 x 75 mm) polystyrene round-bottom tube (e.g. Corning Catalog #352058)	14 mL (17 x 100 mm) polystyrene round-bottom tube (e.g. Corning Catalog #352057)
4	Add Isolation Cocktail to sample.	50 µL/mL of sample	50 µL/mL of sample
	Mix and incubate.	2 - 8°C for 10 minutes	2 - 8°C for 10 minutes
5	Vortex RapidSpheres™.	30 seconds	30 seconds
6	Add RapidSpheres™ to sample.	75 µL/mL of sample	75 µL/mL of sample
	Mix and incubate.	2 - 8°C for 5 minutes	2 - 8°C for 5 minutes
7	Add recommended medium to top up the sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.	Top up to 2.5 mL	<ul style="list-style-type: none"> • Top up to 5 mL for samples ≤ 4 mL • Top up to 10 mL for samples > 4 mL
	Place the tube (without lid) into the magnet and incubate.	<ul style="list-style-type: none"> • For bone marrow, spleen: RT for 5 minutes • For blood: RT for 10 minutes 	<ul style="list-style-type: none"> • For bone marrow, spleen: RT for 5 minutes • For blood: RT for 10 minutes
8	Pick up the magnet, and in one continuous motion invert the magnet and tube,* pouring off the enriched cell suspension into a new tube.	Use new 5 mL tube	Use new 14 mL tube
9	Remove the original tube from the magnet and place the new tube (without lid) into the magnet and incubate for a second separation.	RT for 5 minutes	RT for 5 minutes
10	Pick up the magnet, and in one continuous motion invert the magnet and tube,* pouring off the enriched cell suspension into a new tube.	Isolated cells are ready for use	Isolated cells are ready for use

RT - room temperature (15 - 25°C)

* Leave the magnet and tube inverted for 2 - 3 seconds, then return upright. Do not shake or blot off any drops that may remain hanging from the mouth of the tube.

Table 2. EasySep™ Mouse/Human Chimera Isolation Kit Protocol

STEP	INSTRUCTIONS	EASYSEPT™ MAGNETS		
		 EasyPlate™ (Catalog #18102)	EasyEights™ (Catalog #18103)	
			 5 mL tube	 14 mL tube
1	Prepare sample at the indicated cell concentration within the volume range.	1 x 10 ⁸ cells/mL 0.05 - 0.2 mL	1 x 10 ⁸ cells/mL 0.25 - 2 mL	1 x 10 ⁸ cells/mL 0.5 - 8 mL
2	Add Rat Serum to sample.	50 µL/mL of sample	50 µL/mL of sample	50 µL/mL of sample
3	Add sample to required tube (or plate when using the EasyPlate™ EasySep™ Magnet).	Round bottom, non-tissue culture-treated 96-well plate (e.g. Costar Catalog #3788)	5 mL (12 x 75 mm) polystyrene round-bottom tube (e.g. Corning Catalog #352058)	14 mL (17 x 100 mm) polystyrene round-bottom tube (e.g. Corning Catalog #352057)
4	Add Isolation Cocktail to sample.	50 µL/mL of sample	50 µL/mL of sample	50 µL/mL of sample
	Mix and incubate.	2 - 8°C for 10 minutes	2 - 8°C for 10 minutes	2 - 8°C for 10 minutes
5	Vortex RapidSpheres™.	30 seconds	30 seconds	30 seconds
6	Add RapidSpheres™ to sample.	75 µL/mL of sample	75 µL/mL of sample	75 µL/mL of sample
	Mix and incubate.	2 - 8°C for 5 minutes	2 - 8°C for 5 minutes	2 - 8°C for 5 minutes
7	Add recommended medium to top up sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.	Top up to 0.25 mL	Top up to 2.5 mL	<ul style="list-style-type: none"> • Top up to 5 mL for samples ≤ 4 mL • Top up to 10 mL for samples > 4 mL
	Place the tube or plate (without lid) into the magnet and incubate.	RT for 10 minutes	<ul style="list-style-type: none"> • For bone marrow, spleen: RT for 5 minutes • For blood: RT for 10 minutes 	<ul style="list-style-type: none"> • For bone marrow, spleen: RT for 5 minutes • For blood: RT for 10 minutes
8	Carefully pipette** (do not pour) off the enriched cell suspension into a new tube or plate.	Use new 96-well plate	Use new 5 mL tube	Use new 14 mL tube
9	Remove the tube or plate from the magnet. Place the new tube or plate (without lid) into the magnet and incubate for a second separation.	RT for 5 minutes	RT for 5 minutes	RT for 5 minutes
10	Carefully pipette** (do not pour) off the enriched cell suspension into a new tube or plate.	Isolated cells are ready for use	Isolated cells are ready for use	Isolated cells are ready for use

RT - room temperature (15 - 25°C)

** Collect the entire supernatant, all at once, into a single pipette (e.g. for the EasyEights™ 5 mL tube use a 2 mL serological pipette and for the EasyEights™ 14 mL tube use a 10 mL serological pipette).

Notes and Tips

ASSESSING PURITY


For purity assessment of human cells by flow cytometry use the following fluorochrome-conjugated antibody clone:

- Anti-Human CD45 Antibody, Clone HI30 (Catalog #60018)

For purity assessment of mouse cells by flow cytometry use the following fluorochrome-conjugated antibody clones:

- Anti-Mouse CD45 Antibody, Clone 30-F11 (Catalog #60030), and
- Anti-Mouse TER119 Antibody, Clone TER-119 (Catalog #60033)

Data

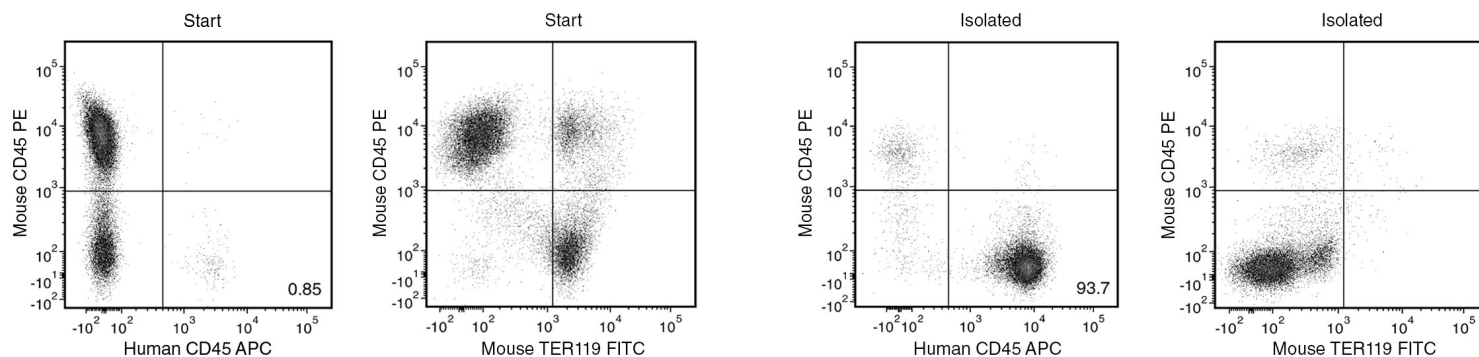

In the above example, human peripheral mononuclear cells (PBMCs) were seeded into mouse bone marrow at a frequency of 0.85%. The human PBMC content of the isolated fraction was 93.7%. The depletion of CD45+ mouse cells is typically 3 - 4.5 log.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, Scientists Helping Scientists, EasyEights, EasyPlate, EasySep, RapidSpheres, and RoboSep are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.