


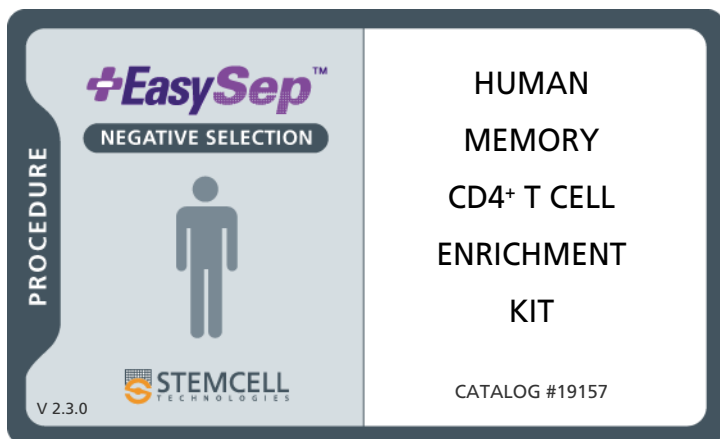

THIS PRODUCT INFORMATION SHEET IS PROVIDED FOR USE WITH ROBOSEP™ (SECTION A), THE PURPLE EASYSEP™ MAGNET (SECTION B) OR "THE BIG EASY" SILVER EASYSEP™ MAGNET (SECTION C).


If using other EasySep™ Magnets, please visit www.stemcell.com to download the magnet-specific Product Information Sheet or contact Technical Support at techsupport@stemcell.com.

A) FULLY AUTOMATED PROTOCOL USING ROBOSEP™ (CATALOG #20000).

This procedure is used for processing **500 µL - 8.5 mL** of sample (up to 4.25×10^8 cells).

1. Prepare cell suspension at a concentration of 5×10^7 cells/mL in RoboSep™ Buffer (Catalog #20104). Cells must be placed in a 14 mL (17 x 100 mm) polystyrene tube to properly fit into the RoboSep™ carousel.
Falcon™ 14 mL Polystyrene Round-Bottom Tubes (BD Biosciences, Catalog #352057) are recommended.
2. Select the appropriate RoboSep™ protocol:
 - Human Memory CD4+ T Cell Negative Selection 19157 - high recovery
 If a modified RoboSep™ protocol is required, please contact STEMCELL Technologies' Technical Support at techsupport@stemcell.com.
3. Load the RoboSep™ carousel as directed by the on-screen prompts. **Vortex the EasySep™ D Magnetic Particles for 30 seconds before loading. Ensure that the particles are in a uniform suspension with no visible aggregates.** When all desired quadrants are loaded, press the green "Run" button. All cell labeling and separation steps will be performed by RoboSep™.
4. When cell separation is complete, remove the enriched cells in the 50 mL tube located to the left of the tip rack. The enriched cells are now ready for use.

MANUAL EASYSEP™ PROTOCOL DIAGRAM

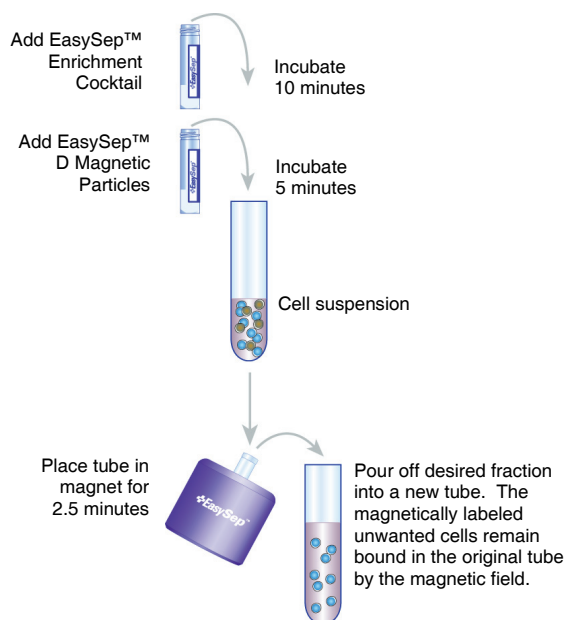

B) MANUAL EASYSEP™ PROTOCOL USING THE PURPLE EASYSEP™ MAGNET (CATALOG #18000).

This procedure is used for processing **250 µL - 2 mL** of sample (up to 1×10^8 cells).

1. Prepare cell suspension at a concentration of 5×10^7 cells/mL in recommended medium (see Notes and Tips, reverse side). Cells must be placed in a 5 mL (12 x 75 mm) polystyrene tube to properly fit into the Purple EasySep™ Magnet.
Falcon™ 5 mL Polystyrene Round-Bottom Tubes (BD Biosciences, Catalog #352058) are recommended.
2. Add the EasySep™ Human Memory CD4+ T Cell Enrichment Cocktail at **50 µL/mL cells** (e.g. for 2 mL of cells, add 100 µL of cocktail). Mix well and incubate at room temperature (15 - 25°C) for **10 minutes**.
3. Vortex the EasySep™ D Magnetic Particles for 30 seconds. Ensure that the particles are in a uniform suspension with no visible aggregates.
4. Add the EasySep™ D Magnetic Particles at **50 µL/mL cells** (e.g. for 2 mL of cells, add 100 µL of magnetic particles). Mix well and incubate at room temperature (15 - 25°C) for **5 minutes**.
5. Bring the cell suspension up to a **total volume of 2.5 mL** by adding recommended medium. Mix the cells in the tube by gently pipetting up and down 2 - 3 times. Place the tube (without cap) into the magnet. Set aside for **2.5 minutes**.
6. Pick up the EasySep™ Magnet, and in one continuous motion invert the magnet and tube, pouring off the desired fraction into a new 5 mL polystyrene tube. The magnetically labeled unwanted cells will remain bound inside the original tube, held by the magnetic field of the EasySep™ Magnet. Leave the magnet and tube in inverted position for 2 - 3 seconds, then return to upright position. *Do not shake or blot off any drops that may remain hanging from the mouth of the tube.* The negatively selected, enriched cells in the new tube are now ready for use.

Additional Notes:

Recommendations for increased recovery: Some of the desired cells may be left behind in the original tube after pouring off the desired fraction. These cells may be recovered by resuspending the magnetically labeled cells in 2.5 mL of recommended medium and performing an additional round of magnetic separation. The supernatant of this separation step can be combined with the primary desired fraction. This will improve recovery, but may decrease purity. Purity assessment prior to pooling is therefore recommended, especially for cell populations with low start percentages.

C) MANUAL EASYSEP™ PROTOCOL USING "THE BIG EASY" SILVER EASYSEP™ MAGNET (CATALOG #18001).

This procedure is used for processing **500 µL - 8.5 mL** of sample (up to 4.25×10^8 cells).

1. Prepare cell suspension at a concentration of 5×10^7 cells/mL in recommended medium (see Notes and Tips, reverse side). Cells must be placed in a 14 mL (17 x 100 mm) polystyrene tube to properly fit into the Silver EasySep™ Magnet.
Falcon™ 14 mL Polystyrene Round-Bottom Tubes (BD Biosciences, Catalog #352057) are recommended.
2. Add the EasySep™ Human Memory CD4+ T Cell Enrichment Cocktail at **50 µL/mL cells** (e.g. for 2 mL of cells, add 100 µL of cocktail). Mix well and incubate at room temperature (15 - 25°C) for **10 minutes**.
3. Vortex the EasySep™ D Magnetic Particles for 30 seconds. Ensure that the particles are in a uniform suspension with no visible aggregates.
4. Add the EasySep™ D Magnetic Particles at **50 µL/mL cells** (e.g. for 2 mL of cells, add 100 µL of magnetic particles). Mix well and incubate at room temperature (15 - 25°C) for **5 minutes**.
5. Bring the cell suspension up to a **total volume of 5 mL** (for $<10^8$ cells) or **10 mL** (for $1 - 4.25 \times 10^8$ cells) by adding recommended medium. Mix the cells in the tube by gently pipetting up and down 2 - 3 times. Place the tube (without cap) into the magnet. Set aside for **2.5 minutes**.
6. Pick up the EasySep™ Magnet, and in one continuous motion invert the magnet and tube, pouring off the desired fraction into a new 14 mL polystyrene tube. The magnetically labeled unwanted cells will remain bound inside the original tube, held by the magnetic field of the EasySep™ Magnet. Leave the magnet and tube in inverted position for 2 - 3 seconds, then return to upright position. *Do not shake or blot off any drops that may remain hanging from the mouth of the tube.* The negatively selected, enriched cells in the new tube are now ready for use.

Additional Notes:

Recommendations for increased recovery: Some of the desired cells may be left behind in the original tube after pouring off the desired fraction. These cells may be recovered by resuspending the magnetically labeled cells in 5 mL (for $<10^8$ cells) or 10 mL (for $1 - 4.25 \times 10^8$ cells) of recommended medium and performing an additional magnetic separation. The supernatant of this separation step can be combined with the primary desired fraction. This will improve recovery, but may decrease purity. Purity assessment prior to pooling is therefore recommended, especially for cell populations with low start percentages.

WWW.STEMCELL.COM

FOR RESEARCH USE ONLY. NOT FOR THERAPEUTIC OR DIAGNOSTIC USE.

Components:

- EasySep™ Human Memory CD4⁺ T Cell Enrichment Cocktail 1.0 mL
- EasySep™ D Magnetic Particles 2 x 1.0 mL



NEGATIVE SELECTION

REQUIRED EQUIPMENT:

EasySep™ Magnet (Catalog #18000), or "The Big Easy" EasySep™ Magnet (Catalog #18001), or RoboSep™ (Catalog #20000).

PRODUCT DESCRIPTION AND APPLICATIONS:

EasySep™ Human Memory CD4⁺ T Cell Enrichment Cocktail and EasySep™ D Magnetic Particles label non-CD4⁺ T cells and naïve CD4⁺ T cells for magnetic separation. These reagents are designed to enrich memory CD4⁺ T cells from fresh or previously frozen peripheral blood mononuclear cells by depletion of non-CD4⁺ T cells and naïve CD4⁺ T cells.

EASYSEPTM LABELING OF HUMAN CELLS:

Unwanted cells are specifically labeled with dextran-coated magnetic particles using bispecific Tetrameric Antibody Complexes (TAC). These complexes recognize both dextran and the cell surface antigen expressed on the unwanted cells (Figure 1). The small size of the magnetic dextran iron particles allows for efficient binding to the TAC-labeled cells. Magnetically labeled cells are then separated from unlabeled target cells using the EasySep™ procedure (reverse side).

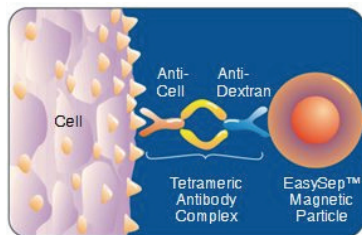


Figure 1.
Schematic Drawing of EasySep™
TAC Magnetic Labeling of Human
Cells.

NOTES AND TIPS:**PREPARING THE CELL SUSPENSION****FROM WHOLE PERIPHERAL BLOOD**

Prepare a mononuclear cell suspension from whole peripheral blood by density gradient centrifugation. **For previously frozen mononuclear cells, we recommend incubating the cells with DNaseI (Catalog #07900)** at a concentration of 100 µg/mL for at least 15 minutes at room temperature (15 - 25°C) prior to labeling and separation. Filter clumpy suspensions through a 30 µm mesh nylon strainer for optimal results.


OPTIMAL CELL NUMBER. The use of fewer than 5×10^7 cells per separation may result in sub-optimal performance.

RECOMMENDED MEDIUM. The recommended medium is RoboSep™ Buffer (Catalog #20104), or Phosphate Buffered Saline (PBS) + 2% FBS (Catalog #07905) with 1 mM EDTA. Medium should be Ca⁺⁺ and Mg⁺⁺ free.

ASSESSING PURITY. Purity of memory CD4⁺ T cells can be measured by flow cytometry after staining with fluorochrome-conjugated anti-CD4 (e.g. APC anti-CD4), anti-CD45RA (e.g. FITC anti-CD45RA, Catalog #10418) and anti-CD45RO (e.g. PE anti-CD45RO, Catalog #10524) antibodies.

TYPICAL EASYSEPTM MEMORY CD4+ T CELL ENRICHMENT PROFILE:

Start: 8% CD4⁺CD45RA⁺CD45RO⁺ Cells Enriched: 93% CD4⁺CD45RA⁺CD45RO⁺ Cells

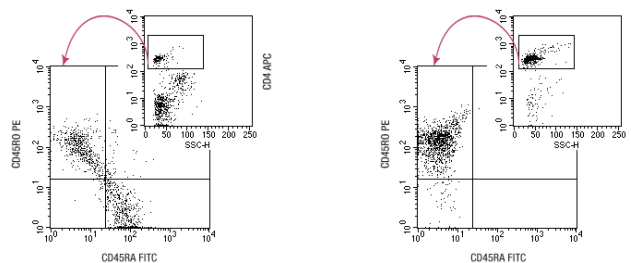

Starting with previously frozen mononuclear cells, the memory CD4⁺ T cell content of the enriched fraction typically ranges from 86 - 98%.

COMPONENT DESCRIPTIONS:**EASYSEPTM HUMAN MEMORY CD4+ T CELL ENRICHMENT COCKTAIL**

CODE #19157C

This cocktail contains a combination of monoclonal antibodies bound in bispecific Tetrameric Antibody Complexes (TAC) which are directed against cell surface antigens on human blood cells (CD8, CD14, CD16, CD19, CD20, CD36, CD45RA, CD56, CD123, TCRγ/δ, glycophorin A) and dextran. The mouse monoclonal antibody subclass is IgG₁. It should be noted that this product is a biological reagent, and as such cannot be completely characterized or quantified. Some variability is unavoidable.

EASYSEPTM D MAGNETIC PARTICLES

CODE #19250

A suspension of magnetic dextran iron particles in TRIS buffer.

STABILITY AND STORAGE:**EASYSEPTM HUMAN MEMORY CD4+ T CELL ENRICHMENT COCKTAIL****EASYSEPTM D MAGNETIC PARTICLES**

Product stable at 2 - 8°C until expiry date as indicated on label. Contents have been sterility tested. Do not freeze this product. This product may be shipped at room temperature (15 - 25°C), and should be refrigerated upon receipt.