



## EasySep™ Human Basophil Enrichment Kit

Negative Selection

Catalog #19069

For processing  $1 \times 10^9$  cells



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Document #29109 | Version 1\_1\_1

## Description

Isolate untouched and highly purified basophils from HetaSep™-processed human peripheral blood by immunomagnetic negative selection.

- Fast, easy-to-use and column-free
- Up to 99% purity
- Isolated cells are untouched

This kit targets non-basophils for removal with antibodies recognizing specific cell surface markers. Unwanted cells are labeled with antibodies and magnetic particles, and separated without columns using an EasySep™ magnet. Desired cells are simply poured off into a new tube. Isolated cells are immediately available for downstream applications such as flow cytometry, culture, or DNA/RNA extraction.

## Component Descriptions

COMPONENT NAME	COMPONENT #	QUANTITY	STORAGE	SHELF LIFE	FORMAT
EasySep™ Human Basophil Enrichment Cocktail	19069C	1 x 1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A combination of monoclonal antibodies in PBS.
EasySep™ Magnetic Nanoparticles	19150.1	3 x 1 mL	Store at 2 - 8°C. Do not freeze.	Stable until expiry date (EXP) on label.	A suspension of magnetic particles in water.

PBS - phosphate-buffered saline

Components may be shipped at room temperature (15 - 25°C) but should be stored as indicated above.

## Sample Preparation

For available fresh and frozen samples, see [www.stemcell.com/primarycells](http://www.stemcell.com/primarycells).

### PERIPHERAL BLOOD

IMPORTANT: Do not use dextran sedimentation to prepare cells.

1. Collect whole blood in a blood collection tube containing anticoagulant.
2. Add 1 part HetaSep™ (Catalog #07906) to 5 parts whole blood and mix well. Use the minimum-sized tube for the total volume of HetaSep™:blood sample. A 14 mL tube is the maximum size recommended for optimal leukocyte recovery.
3. Centrifuge sample at  $110 \times g$  for 6 minutes at room temperature (15 - 25°C) with the brake off.
4. Remove tube from centrifuge and let sit undisturbed until red blood cell (RBC):plasma interface is at approximately 40% of the total volume (a maximum of 10 minutes).
5. Harvest the plasma containing the nucleated cells (everything above the RBC fraction) into a 50 mL tube, and add 4 parts cold recommended medium to 1 part harvested cells/plasma.
6. Centrifuge at  $500 \times g$  for 10 minutes at room temperature with the brake set to low.
7. Discard supernatant and wash pellet to remove excess platelets, centrifuging at  $120 \times g$  for 10 minutes at room temperature with the brake off.
8. Discard supernatant and resuspend cells at  $5 \times 10^7$  cells/mL in recommended medium.



## Recommended Medium

EasySep™ Buffer (Catalog #20144), RoboSep™ Buffer (Catalog #20104), or PBS containing 2% fetal bovine serum (FBS) and 1 mM EDTA. Medium should be free of  $\text{Ca}^{++}$  and  $\text{Mg}^{++}$ .

## Directions for Use – Manual EasySep™ Protocols

See page 1 for Sample Preparation and Recommended Medium. Refer to Table 1 for detailed instructions regarding the EasySep™ procedure for each magnet.

**Table 1. EasySep™ Human Basophil Enrichment Kit Protocol**

		EASYSEP™ MAGNETS	
STEP	INSTRUCTIONS	 <b>EasySep™</b> (Catalog #18000)	 <b>“The Big Easy”</b> (Catalog #18001)
1	Prepare sample at the indicated cell concentration within the volume range.	5 x 10 <sup>7</sup> cells/mL 0.5 - 2 mL	5 x 10 <sup>7</sup> cells/mL 0.5 - 6.5 mL
	Add sample to required tube.	5 mL (12 x 75 mm) polystyrene round-bottom tube (e.g. Corning Catalog #352058)	14 mL (17 x 100 mm) polystyrene round-bottom tube (e.g. Corning Catalog #352057)
2	Add Enrichment Cocktail to sample.	50 µL/mL of sample	50 µL/mL of sample
	Mix and incubate.	RT for 10 minutes	RT for 10 minutes
3	Mix Magnetic Particles. NOTE: Particles should appear evenly dispersed.	Pipette up and down more than 5 times	Pipette up and down more than 5 times
4	Add Magnetic Particles to sample.	100 µL/mL of sample	100 µL/mL of sample
	Mix and incubate.	RT for 10 minutes	RT for 10 minutes
5	Add recommended medium to top up the sample to the indicated volume. Mix by gently pipetting up and down 2 - 3 times.	Top up to 2.5 mL	<ul style="list-style-type: none"> <li>• Top up to 5 mL for samples &lt; 2 mL</li> <li>• Top up to 10 mL for samples ≥ 2 mL</li> </ul>
	Place the tube (without lid) into the magnet and incubate.	RT for 5 minutes	RT for 10 minutes
6	Pick up the magnet, and in one continuous motion invert the magnet and tube,* pouring the enriched cell suspension into a new tube.	Use a new 5 mL tube	Use a new 14 mL tube
7	Remove the tube from the magnet and place the new tube (without lid) into the magnet and incubate.	RT for 5 minutes	RT for 10 minutes
8	Repeat steps as indicated.	Steps 6 and 7 (total of 3 x 5-minute separations)	---
9	Pick up the magnet, and in one continuous motion invert the magnet and tube,* pouring the enriched cell suspension into a new tube.	Isolated cells are ready for use	Isolated cells are ready for use


RT - room temperature (15 - 25°C)

\* Leave the magnet and tube inverted for 2 - 3 seconds, then return upright. Do not shake or blot off any drops that may remain hanging from the mouth of the tube.

## Directions for Use – Fully Automated RoboSep™ Protocol

See page 1 for Sample Preparation and Recommended Medium. Refer to Table 2 for detailed instructions regarding the RoboSep™ procedure.

**Table 2. RoboSep™ Human Basophil Enrichment Kit Protocol**

STEP	INSTRUCTIONS	RoboSep™ (Catalog #20000 and #21000)	
1	Prepare sample at the indicated cell concentration within the volume range.	5 x 10 <sup>7</sup> cells/mL 0.5 - 6.5 mL	
	Add sample to required tube.	14 mL (17 x 100 mm) polystyrene round-bottom tube (e.g. Corning Catalog #352057)	
2	Select protocol.	Human Basophil Negative Selection 19069-high purity	
3	Mix Magnetic Particles. NOTE: Particles should appear evenly dispersed.	Pipette up and down more than 5 times	
4	Load the carousel.	Follow on-screen prompts	
	Start the protocol.	Press the green "Run" button	
5	Unload the carousel when the run is complete.	Isolated cells are ready for use	

## Notes and Tips


### ASSESSING PURITY

For purity assessment of human basophils (CD123+IgE+ or CD123+CD203c<sup>low</sup>) by flow cytometry use the following fluorochrome-conjugated antibody clones:

- Anti-Human CD123 (IL-3Ra) Antibody, Clone 6H6 (Catalog #60110), and
- Anti-human IgE antibody, or
- Anti-human CD203c antibody

Alternatively, purity may be assessed by performing a cytospin on the enriched cells followed by Wright's or May-Grünwald staining (e.g. Sigma-Aldrich Catalog #W0625 or #205435, respectively).

## Data

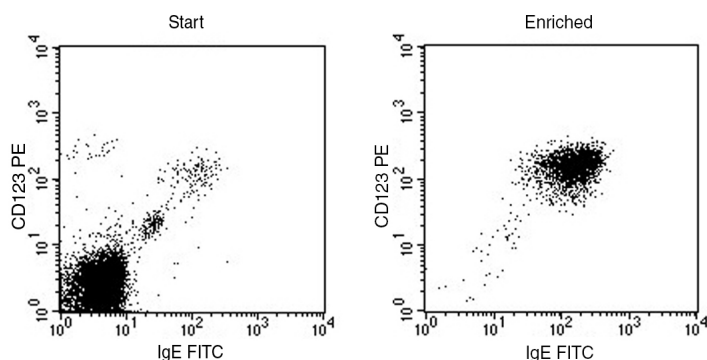

Starting with freshly prepared nucleated cells, the basophil content (CD123+IgE+) of the enriched fraction typically ranges from 92 - 99% (gated on CD45+). In the above example, the purities of the start and final enriched fractions are 1.5% and 98%, respectively.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists, EasySep, HetaSep, and RoboSep are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.