

Small Molecules

1-Oleoyl Lysophosphatidic Acid

Lysophosphatidic acid (LPA)1 and LPA2 agonist

Catalog # 72692
72694

1 mg
5 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

1-Oleoyl Lysophosphatidic Acid is a species of lysophosphatidic acid (LPA) containing oleic acid at the *sn*-1 position. LPA binds to four different G-protein-linked receptors (Chun et al.) to mediate a variety of biological responses including cell proliferation, smooth muscle contraction, platelet aggregation, neurite retraction, and cell motility (Moolenaar). 1-Oleoyl Lysophosphatidic Acid is the most potent of the LPA analogs for calcium mobilization in A431 cells (Jalink et al.) and for growth stimulation of a variety of cell lines (van Corven et al.). This product is supplied as the sodium salt of the molecule.

Molecular Name:	1-Oleoyl Lysophosphatidic Acid (Sodium Salt)
Alternative Names:	Oleoyl- <i>sn</i> -3-glycerophosphate
CAS Number:	22556-62-3
Chemical Formula:	C ₂₁ H ₄₀ O ₇ P · Na
Molecular Weight:	458.5 g/mol
Purity:	≥ 95%
Chemical Name:	sodium;[2-hydroxy-3-[(Z)-octadec-9-enoyl]oxypropyl] hydrogen phosphate
Structure:	

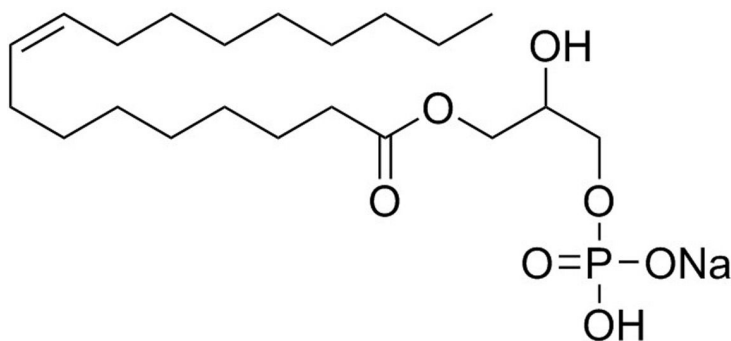

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com .
Solubility:	· PBS (pH 7.2) ≤ 18 mM For example, to prepare a 10 mM stock solution in PBS, resuspend 1 mg in 218 µL of PBS (pH 7.2).

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in PBS at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use.

Published Applications

DIFFERENTIATION

- Stimulates neuronal differentiation in cultured mouse or rat neural progenitor cells (Cui & Qiao; Fukushima et al.; Spohr et al.).
- Inhibits human embryonic stem (ES) cell-derived neural stem cells (NSCs) from forming neurospheres and differentiating into neurons in vitro (Dottori et al.).
- Stimulates differentiation of human adipose tissue-derived mesenchymal stem cells to myofibroblast-like cells in vitro (Jeon et al.).

References

- Chun J et al. (2002) International Union of Pharmacology. XXXIV. Lysophospholipid receptor nomenclature. *Pharmacol Rev* 54(2): 265–9.
- Cui HL & Qiao JT. (2006) Promotive action of lysophosphatidic acid on proliferation of rat embryonic neural stem cells and their differentiation to cholinergic neurons in vitro. *Sheng Li Xue Bao* 58(6): 547–55.
- Dottori M et al. (2008) Lysophosphatidic acid inhibits neuronal differentiation of neural stem/progenitor cells derived from human embryonic stem cells. *Stem Cells* 26(5): 1146–54.
- Fukushima N et al. (2007) Lysophosphatidic acid stimulates neuronal differentiation of cortical neuroblasts through the LPA1-G(i/o) pathway. *Neurochem Int* 50(2): 302–7.
- Jalink K et al. (1995) Lysophosphatidic acid-induced Ca²⁺ mobilization in human A431 cells: structure-activity analysis. *Biochem J* 307 (Pt 2): 609–16.
- Jeon ES et al. (2008) Cancer-derived lysophosphatidic acid stimulates differentiation of human mesenchymal stem cells to myofibroblast-like cells. *Stem Cells* 26(3): 789–97.
- Moolenaar WH. (1994) LPA: a novel lipid mediator with diverse biological actions. *Trends Cell Biol* 4(6): 213–9.
- Spohr TCS et al. (2008) Lysophosphatidic acid receptor-dependent secondary effects via astrocytes promote neuronal differentiation. *J Biol Chem* 283(12): 7470–9.
- Van Corven EJ et al. (1992) Mitogenic action of lysophosphatidic acid and phosphatidic acid on fibroblasts. Dependence on acyl-chain length and inhibition by suramin. *Biochem J* 281 (Pt 1): 163–9.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.