

Small Molecules

GW2580

CSF-1 pathway inhibitor; Inhibits cFMS

Catalog # 72472
72474

5 mg
25 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

GW2580 is a selective inhibitor of cFMS kinase ($IC_{50} = 0.03 \mu M$; Conway et al. 2008), blocking its ability to autophosphorylate colony-stimulating factor (CSF-1 or M-CSF), which promotes the survival, proliferation, and differentiation of macrophages.

Molecular Name:	GW2580
Alternative Names:	Not applicable
CAS Number:	870483-87-7
Chemical Formula:	$C_{20}H_{22}N_4O_3$
Molecular Weight:	366.4 g/mol
Purity:	$\geq 98\%$
Chemical Name:	5-[[[3-methoxy-4-[(4-methoxyphenyl)methoxy]phenyl]methyl]-2,4-pyrimidinediamine
Structure:	

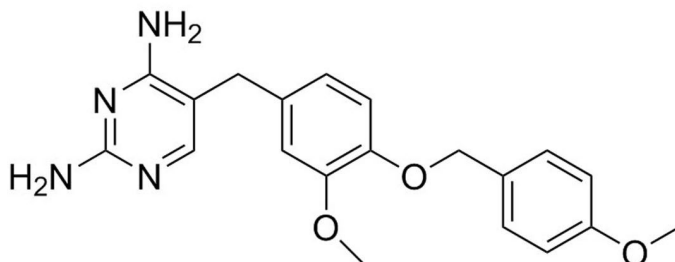

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at $-20^{\circ}C$ as supplied. Protect from prolonged exposure to light. Stable as supplied for 12 months from date of receipt.
Solubility:	· DMSO ≤ 25 mM For example, to prepare a 1 mM stock solution in DMSO, resuspend 1 mg in 2.73 mL of fresh DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at $-20^{\circ}C$. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

DIFFERENTIATION

- Demonstrates the importance of CSF-1 in promoting myeloid lineage bias in mouse hematopoietic stem cells (Mossadegh-Keller et al.).
- Demonstrates the importance of CSF-1 in expansion of mouse and human macrophage colonies and monocytes (He et al.; Clanchy & Hamilton; Conway et al. 2008; Conway et al. 2005).
- Inhibits bone degradation in cultures of human osteoclasts, rat calvaria, and rat fetal long bone (Conway et al. 2005), and in mouse models of arthritis (Conway et al. 2008).

References

- Clanchy FIL & Hamilton JA. (2012) HUVEC co-culture and haematopoietic growth factors modulate human proliferative monocyte activity. *Cytokine* 59(1): 31–4.
- Conway JG et al. (2008) Effects of the cFMS kinase inhibitor 5-(3-methoxy-4-((4-methoxybenzyl)oxy)benzyl)pyrimidine-2,4-diamine (GW2580) in normal and arthritic rats. *J Pharmacol Exp Ther* 326(1): 41–50.
- Conway JG et al. (2005) Inhibition of colony-stimulating-factor-1 signaling in vivo with the orally bioavailable cFMS kinase inhibitor GW2580. *Proc Natl Acad Sci USA* 102(44): 16078–83.
- He H et al. (2012) Endothelial cells provide an instructive niche for the differentiation and functional polarization of M2-like macrophages. *Blood* 120(15): 3152–62.
- Mossadegh-Keller N et al. (2013) M-CSF instructs myeloid lineage fate in single haematopoietic stem cells. *Nature* 497(7448): 239–43.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.