

## Small Molecules

KY02111

WNT pathway inhibitor

Catalog # 72582

5 mg



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

KY02111 promotes the differentiation of pluripotent stem cells (PSCs) to cardiomyocytes by inhibiting canonical WNT signaling in a manner distinct from that of other WNT inhibitors (Minami et al.).

Molecular Name:	KY02111
Alternative Names:	Not applicable
CAS Number:	1118807-13-8
Chemical Formula:	C <sub>18</sub> H <sub>17</sub> ClN <sub>2</sub> O <sub>3</sub> S
Molecular Weight:	376.9 g/mol
Purity:	≥ 95%
Chemical Name:	N-(6-chloro-1,3-benzothiazol-2-yl)-3-(3,4-dimethoxyphenyl)propanamide
Structure:	

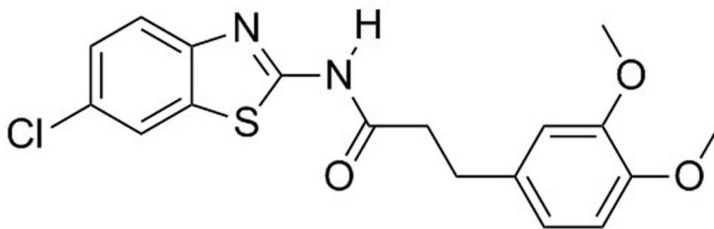

## Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact <a href="mailto:techsupport@stemcell.com">techsupport@stemcell.com</a> .
Solubility:	· DMSO ≤ 25 mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 265 µL of fresh DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

## Published Applications

### DIFFERENTIATION

- Promotes cardiomyocyte differentiation of human and mouse PSCs in combination with BIO, CHIR99021, and XAV939 (Minami et al.).

## References

Minami I et al. (2012) A small molecule that promotes cardiac differentiation of human pluripotent stem cells under defined, cytokine- and xeno-free conditions. Cell Rep 2(5): 1448–60.

## Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at [www.stemcell.com/smallmolecules](http://www.stemcell.com/smallmolecules) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

**This product is hazardous. Please refer to the Safety Data Sheet (SDS).**

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.