

Small Molecules

Rosiglitazone

PPAR γ activator

Catalog # 72622
72624

10 mg
50 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Rosiglitazone is a potent and selective ligand for peroxisome proliferator-activated receptor gamma (PPAR γ). It binds to the PPAR γ ligand-binding domain with a K_d of 43 nM (Lehmann et al.). It activates luciferase-based expression constructs PPAR γ 1 and PPAR γ 2 with EC_{50} values of approximately 30 nM and 100 nM, respectively (Lehmann et al.).

Molecular Name:	Rosiglitazone
Alternative Names:	Avandia; BRL 49653
CAS Number:	122320-73-4
Chemical Formula:	$C_{18}H_{19}N_3O_3S$
Molecular Weight:	357.4 g/mol
Purity:	$\geq 98\%$
Chemical Name:	5-[[4-[2-(methyl-2-pyridinylamino)ethoxy]phenyl]methyl]-2,4-thiazolidinedione
Structure:	

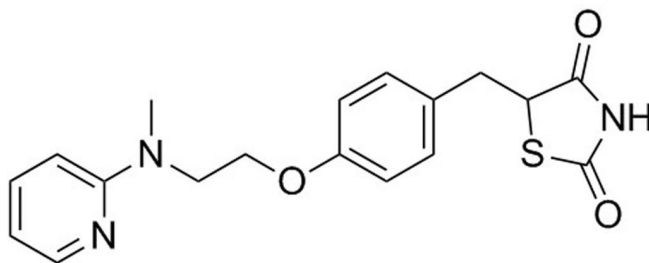

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com .
Solubility:	· Absolute ethanol ≤ 2.7 mM · DMSO ≤ 95 mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 10 mg in 2.8 mL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C . Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

MAINTENANCE AND SELF-RENEWAL

- Stimulates embryonic mouse neural stem cell proliferation and inhibits neuronal differentiation (Wada et al.).

DIFFERENTIATION

- Induces adipocyte differentiation in C3H10T1/2 stem cells (Lehmann et al.).
- Promotes endothelial differentiation and inhibits smooth muscle differentiation in angiogenic progenitor cells (Wang et al.).
- Stimulates adipocyte differentiation and decreases osteogenesis in mouse and human bone marrow-derived mesenchymal stem cells (Ali et al.; Benvenuti et al.; Lecka-Czernik et al.; Sorocéanu et al.).
- Enhances osteoclast formation in mouse bone marrow cells stimulated with RANKL and M-CSF (Wu et al.).

References

- Ali AA et al. (2005) Rosiglitazone causes bone loss in mice by suppressing osteoblast differentiation and bone formation. *Endocrinology* 146(3): 1226–35.
- Benvenuti S et al. (2007) Rosiglitazone stimulates adipogenesis and decreases osteoblastogenesis in human mesenchymal stem cells. *J Endocrinol Invest* 30(9): RC26–30.
- Lecka-Czernik B et al. (2002) Divergent effects of selective peroxisome proliferator-activated receptor-gamma 2 ligands on adipocyte versus osteoblast differentiation. *Endocrinology* 143(6): 2376–84.
- Lehmann JM et al. (1995) An Antidiabetic Thiazolidinedione Is a High Affinity Ligand for Peroxisome Proliferator-activated Receptor (PPAR). *J Biol Chem* 270(22): 12953–6.
- Sorocéanu MA et al. (2004) Rosiglitazone impacts negatively on bone by promoting osteoblast/osteocyte apoptosis. *J Endocrinol* 183(1): 203–16.
- Wada K et al. (2006) Peroxisome proliferator-activated receptor gamma-mediated regulation of neural stem cell proliferation and differentiation. *J Biol Chem* 281(18): 12673–81.
- Wang C-H et al. (2004) Rosiglitazone facilitates angiogenic progenitor cell differentiation toward endothelial lineage: a new paradigm in glitazone pleiotropy. *Circulation* 109(11): 1392–400.
- Wu H et al. (2013) Regulation of selective PPAR γ modulators in the differentiation of osteoclasts. *J Cell Biochem* 114(9): 1969–77.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.