

Small Molecules

Pluripotin

PI3K and MEK/ERK pathway inhibitor;
Inhibits RasGAP and ERK1

Catalog # 72812
72814

1 mg
10 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM


FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Pluripotin inhibits both Ras GTPase-activating protein (RasGAP, $K_d = 98$ nM) and extracellular signal-regulated kinase 1 (ERK1, $K_d = 212$ nM; Chen et al.). Inhibition of RasGAP increases Ras signaling via the PI3-kinase pathway, which promotes self-renewal, whereas inhibition of ERK1 blocks differentiation (Chen et al.).

Molecular Name:	Pluripotin
Alternative Names:	SC-1
CAS Number:	839707-37-8
Chemical Formula:	$C_{27}H_{25}F_3N_8O_2$
Molecular Weight:	550.5 g/mol
Purity:	$\geq 98\%$
Chemical Name:	N-(3-(7-(1,3-dimethyl-1H-pyrazol-5-ylamino)-1-methyl-2-oxo-1,2-dihydropyrimido[4,5-d]pyrimidin-3(4H)-yl)-4-methylphenyl)-3-(trifluoromethyl)benzamide

Structure:

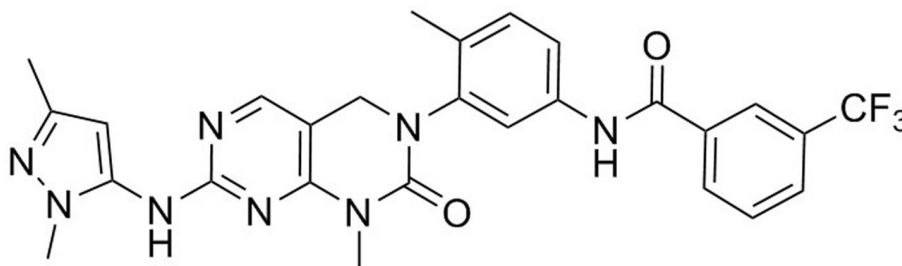

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. Stable as supplied for 12 months from date of receipt.
Solubility:	\cdot DMSO ≤ 45 mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 182 μL of fresh DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C . Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

MAINTENANCE AND SELF-RENEWAL

- Promotes self-renewal and maintains mouse embryonic stem (ES) cells in an undifferentiated, pluripotent state in the absence of feeder cells, serum, and leukemia inhibitory factor (LIF) for at least 10 passages (Chen et al.; Xiong et al.).
- When combined with LIF (Catalog #78056), improves efficiency of ES cell derivation from refractory mouse strains (Pieters et al.; Yang et al.).

References

- Chen S et al. (2006) Self-renewal of embryonic stem cells by a small molecule. *Proc Natl Acad Sci USA* 103(46): 17266–71.
- Pieters T et al. (2012) Efficient and user-friendly pluripotin-based derivation of mouse embryonic stem cells. *Stem Cell Rev* 8(3): 768–78.
- Xiong W et al. (2009) The use of SC1 (pluripotin) to support mESC self-renewal in the absence of LIF. *J Vis Exp* (33): e1550.
- Yang W et al. (2009) Pluripotin combined with leukemia inhibitory factor greatly promotes the derivation of embryonic stem cell lines from refractory strains. *Stem Cells* 27(2): 383–9.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.