

Small Molecules

Pyrintegrin

Integrin and tyrosine kinase activator;
Activates Integrin, FGFR, IGFR, and
HER2

Catalog # 72842
72844

1 mg
10 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM


FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Pyrintegrin is a 2,4-disubstituted pyrimidine that induces the activation of $\beta 1$ integrin and multiple growth factor receptors, including FGFR1, IGFR1, EGFR1, and HER2 (Xu et al.).

Molecular Name:	Pyrintegrin
Alternative Names:	Not applicable
CAS Number:	1228445-38-2
Chemical Formula:	C ₂₃ H ₂₅ N ₅ O ₃ S
Molecular Weight:	451.5 g/mol
Purity:	≥ 95%
Chemical Name:	N-(cyclopropylmethyl)-4-[[4-(3,4-dihydro-6-hydroxy-1(2H)-quinoliny)-2-pyrimidinyl]amino]-benzenesulfonamide

Structure:

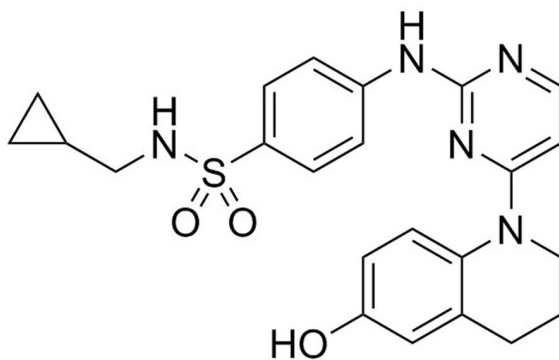

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. Stable as supplied for 12 months from date of receipt.
Solubility:	· Absolute ethanol ≤ 2.2 mM · DMSO ≤ 65 mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 221 μ L of fresh DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

MAINTENANCE AND SELF-RENEWAL

- Enhances integrin-dependent attachment and survival of human embryonic stem cells following trypsin-mediated single-cell dissociation (Xu et al.).

References

Xu Y et al. (2010) Revealing a core signaling regulatory mechanism for pluripotent stem cell survival and self-renewal by small molecules. Proc Natl Acad Sci USA 107(18): 8129–34.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, visit www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.