

## Small Molecules

### Zebularine

Epigenetic modifier; Inhibits DNA methyltransferases (DNMT)

Catalog # 72902

10 mg



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

Zebularine, a cytidine analog, is a DNA methylation inhibitor that acts by forming a covalent complex with DNA methyltransferases (DNMTs; Zhou et al.). By stabilizing the binding of DNMTs to DNA, zebularine hinders methylation and decreases dissociation, trapping DNMT and preventing its turnover (Champion et al.). Zebularine is also a potent inhibitor of cytidine deaminase (Laliberté et al.).

|                    |                                                                             |
|--------------------|-----------------------------------------------------------------------------|
| Molecular Name:    | Zebularine                                                                  |
| Alternative Names: | NSC 309132                                                                  |
| CAS Number:        | 3690-10-6                                                                   |
| Chemical Formula:  | C <sub>9</sub> H <sub>12</sub> N <sub>2</sub> O <sub>5</sub>                |
| Molecular Weight:  | 228.2 g/mol                                                                 |
| Purity:            | ≥ 98%                                                                       |
| Chemical Name:     | 1-[(2R,3R,4S,5R)-3,4-dihydroxy-5-(hydroxymethyl)oxolan-2-yl]pyrimidin-2-one |
| Structure:         |                                                                             |

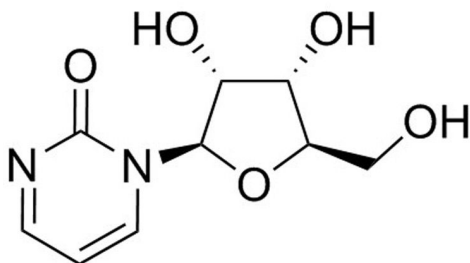

## Properties

|                      |                                                                                                                                                                                                                                    |
|----------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Physical Appearance: | A crystalline solid                                                                                                                                                                                                                |
| Storage:             | Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact <a href="mailto:techsupport@stemcell.com">techsupport@stemcell.com</a> .                                    |
| Solubility:          | <ul style="list-style-type: none"><li>· PBS (pH 7.2) ≤ 40 mM</li><li>· DMSO ≤ 60 mM</li><li>· Absolute ethanol ≤ 1 mM</li></ul> For example, to prepare a 10 mM stock solution in PBS, resuspend 10 mg in 4.38 mL of PBS (pH 7.2). |

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

## Published Applications

### DIFFERENTIATION

- Induces cardiomyocyte differentiation from rat bone marrow mesenchymal stem cells (MSCs; Naeem et al.).

### CANCER RESEARCH

- Reduces proliferation of the cancer cell lines TK6, Jurkat, KG-1, and HCT116 (Stresemann et al.).
- Inhibits cell proliferation in T24 bladder carcinoma cells (Ben-Kasus et al.).
- In combination with the histone deacetylase inhibitor SAHA (Catalog #73902), reduces cell proliferation and increases apoptosis in pancreatic cancer cell lines (Neureiter et al.).
- Inhibits cell proliferation and induces apoptosis in human acute myeloid leukemia cells in vitro (Scott et al.).
- Increases the proportion of cells with cancer stem cell properties in subpopulations of human cancer cells and liver cancer cell lines (Marquardt et al.).

## References

- Ben-Kasus T et al. (2005) Metabolic activation of zebularine, a novel DNA methylation inhibitor, in human bladder carcinoma cells. *Biochem Pharmacol* 70(1): 121–33.
- Champion C et al. (2010) Mechanistic insights on the inhibition of c5 DNA methyltransferases by zebularine. *PLoS One* 5(8): e12388.
- Laliberté J et al. (1992) Potent inhibitors for the deamination of cytosine arabinoside and 5-aza-2'-deoxycytidine by human cytidine deaminase. *Cancer Chemother Pharmacol* 30(1): 7–11.
- Marquardt JU et al. (2011) Human hepatic cancer stem cells are characterized by common stemness traits and diverse oncogenic pathways. *Hepatology* 54(3): 1031–42.
- Naeem N et al. (2013) DNA methylation inhibitors, 5-azacytidine and zebularine potentiate the transdifferentiation of rat bone marrow mesenchymal stem cells into cardiomyocytes. *Cardiovasc Ther* 31(4): 201–9.
- Neureiter D et al. (2007) Apoptosis, proliferation and differentiation patterns are influenced by Zebularine and SAHA in pancreatic cancer models. *Scand J Gastroenterol* 42(1): 103–16.
- Scott SA et al. (2007) Zebularine inhibits human acute myeloid leukemia cell growth in vitro in association with p15INK4B demethylation and reexpression. *Exp Hematol* 35(2): 263–73.
- Stresemann C et al. (2006) Functional diversity of DNA methyltransferase inhibitors in human cancer cell lines. *Cancer Res* 66(5): 2794–800.
- Zhou L et al. (2002) Zebularine: a novel DNA methylation inhibitor that forms a covalent complex with DNA methyltransferases. *J Mol Biol* 321(4): 591–9.

## Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at [www.stemcell.com/smallmolecules](http://www.stemcell.com/smallmolecules) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.