

Small Molecules

Brefeldin A

Protein trafficking inhibitor; Inhibits Sec7-containing guanine-exchange factor (GEF)

Catalog # 73012
73014

10 mg
25 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Brefeldin A is a fungal lactone antibiotic produced by many species, including *Eupenicillium brefeldianum* (Klausner et al.). It reversibly interferes with protein trafficking and secretion mediated by the Golgi apparatus and endoplasmic reticulum by indirect inhibition of ADP-ribosylation factor (ARF; Klausner et al.; Helms & Rothman; Robinson et al.; Morinaga et al.; Moss & Vaughan; Nebenführ et al.; Ktistakis et al.). Brefeldin A binds to Sec7-containing guanine-exchange factor (GEF) at the ARF-GDP-Sec7 interface, preventing the conformational change required to release GDP and activate ARF (Mossessova et al.).

Molecular Name: Brefeldin A

Alternative Names: Ascotoxin, BFA, Cyanein, Decumbin, Nectrolide, NSC 56310, NSC 89671, NSC 107456, NSC 244390, Synergisidin

CAS Number: 20350-15-6


Chemical Formula: $C_{16}H_{24}O_4$

Molecular Weight: 280.4 g/mol

Purity: $\geq 98\%$

Chemical Name: (1R,2E,6S,10E,11aS,13S,14aR)-1,6,7,8,9,11a,12,13,14,14a-decahydro-1,13-dihydroxy-6-methyl-4H-cyclopent[f]oxacyclotridecin-4-one

Structure:

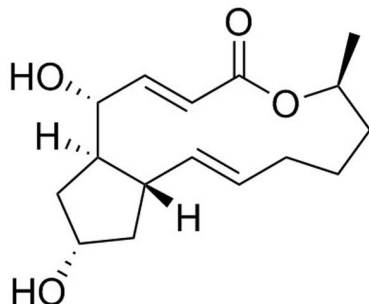

Properties

Physical Appearance: A crystalline solid

Storage: Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com.

Solubility:

- DMSO ≤ 1.5 mM
- Absolute ethanol ≤ 3 mM

For example, to prepare a 1 mM stock solution in DMSO, resuspend 10 mg in 35.7 mL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C . Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

CELL LINE DEVELOPMENT

- Improves clustered regularly interspaced palindromic repeats (CRISPR)-mediated homology-directed repair (HDR) in mouse embryonic stem cells (Yu et al.).

CANCER RESEARCH

- Induces apoptosis in human leukemia (HL60, K562) and colon carcinoma (HT-29) cell lines (Shao et al.).
- Reduces survival, induces apoptosis and inhibits clonogenic activity of Colo 205 colorectal cancer stem cell line (Tseng et al.).

References

- Helms JB & Rothman JE. (1992) Inhibition by brefeldin A of a Golgi membrane enzyme that catalyses exchange of guanine nucleotide bound to ARF. *Nature* 360(6402): 352–4.
- Klausner RD et al. (1992) Brefeldin A: insights into the control of membrane traffic and organelle structure. *J Cell Biol* 116(5): 1071–80.
- Ktistakis NT et al. (1992) Action of brefeldin A blocked by activation of a pertussis-toxin-sensitive G protein. *Nature* 356(6367): 344–6.
- Morinaga N et al. (1999) Brefeldin A inhibited activity of the sec7 domain of p200, a mammalian guanine nucleotide-exchange protein for ADP-ribosylation factors. *J Biol Chem* 274(25): 17417–23.
- Moss J & Vaughan M. (1999) Activation of toxin ADP-ribosyltransferases by eukaryotic ADP-ribosylation factors. *Mol Cell Biochem* 193(1-2): 153–7.
- Mossessova E et al. (2003) Crystal structure of ARF1*Sec7 complexed with Brefeldin A and its implications for the guanine nucleotide exchange mechanism. *Mol Cell* 12(6): 1403–11.
- Nebenführ A et al. (2002) Brefeldin A: deciphering an enigmatic inhibitor of secretion. *Plant Physiol* 130(3): 1102–8.
- Robinson DG et al. (2008) The Endosomal System of Plants: Charting New and Familiar Territories. *PLANT Physiol* 147(4): 1482–1492.
- Shao RG et al. (1996) Brefeldin A is a potent inducer of apoptosis in human cancer cells independently of p53. *Exp Cell Res* 227(2): 190–6.
- Tseng C-N et al. (2013) Brefeldin a effectively inhibits cancer stem cell-like properties and MMP-9 activity in human colorectal cancer Colo 205 cells. *Molecules* 18(9): 10242–53.
- Yu C et al. (2015) Small Molecules Enhance CRISPR Genome Editing in Pluripotent Stem Cells. *Cell Stem Cell* 16(2): 142–147.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.