

Small Molecules

Everolimus

mTOR pathway inhibitor; Inhibits FKBP-12

Catalog # 73122
73124

10 mg
25 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM


FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Everolimus is an inhibitor of mammalian target of rapamycin (mTOR). It is a hydroxyethyl ether substituted derivative of Rapamycin with improved pharmacokinetic and pharmacodynamic properties. It inhibits both mTORC1 and mTORC2 complexes by binding to FK506 Binding Protein (FKBP-12), which then binds to mTOR, leading to complex destabilization and blocked kinase function (Zeng et al.; Lebwohl et al.; Sedrani et al.; Huang & Houghton).

Molecular Name: Everolimus
Alternative Names: Afinitor™, Certican, NVP-RAD001, RAD001, SDZ-RAD, Xience, Zortress™
CAS Number: 159351-69-6
Chemical Formula: C₅₃H₈₃NO₁₄
Molecular Weight: 958.2 g/mol
Purity: ≥ 98%
Chemical Name: (1R,9S,12S,15R,16E,18R,19R,21R,23S,24E,26E,28E,30S,32S,35R)-1,18-Dihydroxy-12-((1R)-2-((1S,3R,4R)-4-(2-hydroxyethoxy)-3-methoxycyclohexyl)-1-methylethyl)-19,30-dimethoxy-15,17,21,23,29,35-hexamethyl-11,36-dioxo-4-azatricyclo(30.3.1.0(sup 4,9))hexatriacont

Structure:

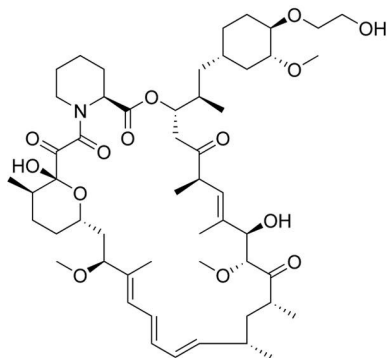

Properties

Physical Appearance: A crystalline solid
Storage: Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com.
Solubility:
· DMSO ≤ 10 mM
· Absolute ethanol ≤ 10 mM
For example, to prepare a 5 mM stock solution in DMSO, resuspend 10 mg in 2.09 mL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

IMMUNOLOGY

- Acts as an immunosuppressive agent in the context of organ transplantation (Lebwohl et al.; Wullschlegler et al.).

CANCER RESEARCH

- Inhibits cell proliferation, metabolism, and angiogenesis in a variety of cancers using in vitro and in vivo models (Zhu et al.; Lane et al.; O'Reilly et al.; Lebwohl et al.).

References

Huang S & Houghton PJ. (2003) Targeting mTOR signaling for cancer therapy. *Curr Opin Pharmacol* 3(4): 371–7.

Lane HA et al. (2009) mTOR inhibitor RAD001 (everolimus) has antiangiogenic/vascular properties distinct from a VEGFR tyrosine kinase inhibitor. *Clin Cancer Res* 15(5): 1612–22.

Lebwohl D et al. (2013) Development of everolimus, a novel oral mTOR inhibitor, across a spectrum of diseases. *Ann N Y Acad Sci* 1291: 14–32.

O'Reilly T et al. (2011) Evaluation of the mTOR inhibitor, everolimus, in combination with cytotoxic antitumor agents using human tumor models in vitro and in vivo. *Anticancer Drugs* 22(1): 58–78.

Sedrani R et al. (1998) Chemical modification of rapamycin: the discovery of SDZ RAD. *Transplant Proc* 30(5): 2192–4.

Wullschlegler S et al. (2006) TOR signaling in growth and metabolism. *Cell* 124(3): 471–84.

Zeng Z et al. (2007) Rapamycin derivatives reduce mTORC2 signaling and inhibit AKT activation in AML. *Blood* 109(8): 3509–12.

Zhu Y et al. (2012) Antitumor effect of the mTOR inhibitor everolimus in combination with trastuzumab on human breast cancer stem cells in vitro and in vivo. *Tumour Biol* 33(5): 1349–62.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.