

Small Molecules

GSA 10

Hedgehog pathway activator;
Activates Smoothened (SMO)

Catalog # 73172
73174

10 mg
50 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

GSA 10 is an agonist of smoothened (SMO), a cell surface receptor and mediator of the hedgehog signaling pathway. It is a quinolinecarboxamide derivative that binds in a distinct binding pocket from cyclopamine. GSA10 is active at an EC_{50} of 1.2 μ M in an alkaline phosphatase-based mouse mesenchymal cell line (C3H10T1/2) differentiation assay (Gorojankina et al.).

Molecular Name:	GSA 10
Alternative Names:	Not applicable
CAS Number:	300833-95-8
Chemical Formula:	$C_{26}H_{30}N_2O_5$
Molecular Weight:	450.5 g/mol
Purity:	≥ 95%
Chemical Name:	4-[[[(1-hexyl-1,2-dihydro-4-hydroxy-2-oxo-3-quinolinyl)carbonyl]amino]-benzoic acid, propyl ester
Structure:	

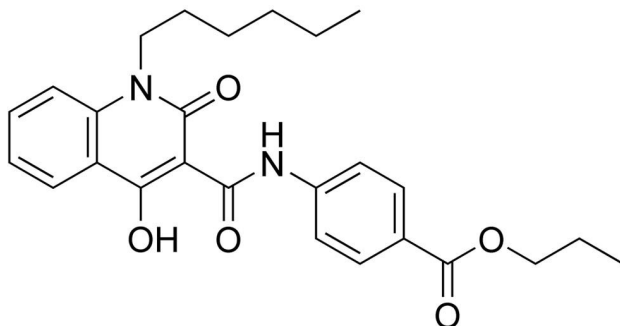

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com .
Solubility:	· DMF ≤ 2 mM For example, to prepare a 1 mM stock solution in DMF, resuspend 10 mg in 22.2 mL of DMF.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMF at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMF concentration above 0.1% due to potential cell toxicity.

Published Applications

DIFFERENTIATION

- Promotes differentiation of multipotent mesenchymal progenitor cells into osteoblasts (Gorojankina et al.).

References

Gorojankina T et al. (2013) Discovery, molecular and pharmacological characterization of GSA-10, a novel small-molecule positive modulator of Smoothened. Mol Pharmacol 83(5): 1020–9.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.