

Small Molecules

(S)-MG132

Ubiquitin-proteasome and NF- κ B pathway inhibitor; Inhibits proteasome and stabilizes I κ B

Catalog # 73262
73264

1 mg
10 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

(S)-MG132 is a reversible, cell-permeable inhibitor of proteasome activity (IC_{50} = 100 nM; Kisselev & Goldberg) and calpain (IC_{50} = 1.2 μ M; Tsubuki et al.). The ubiquitin-proteasome pathway selectively degrades intracellular proteins, thereby clearing damaged or misfolded proteins, and regulating the availability of key proteins involved in the control of inflammatory processes and cell cycle regulation. (S)-MG132 suppresses NF- κ B activation by preventing I κ B degradation (IC_{50} = 3 μ M; Arlt et al., Palombella et al., Ortiz-Lazareno et al.).

Molecular Name:	(S)-MG132
Alternative Names:	Z-Leu-Leu-Leu-CHO; Calpain inhibitor IV-2
CAS Number:	133407-82-6
Chemical Formula:	C ₂₆ H ₄₁ N ₃ O ₅
Molecular Weight:	475.6 g/mol
Purity:	≥ 98%
Chemical Name:	N-[(phenylmethoxy)carbonyl]-L-leucyl-N-[(1S)-1-formyl-3-methylbutyl]-L-leucinamide
Structure:	

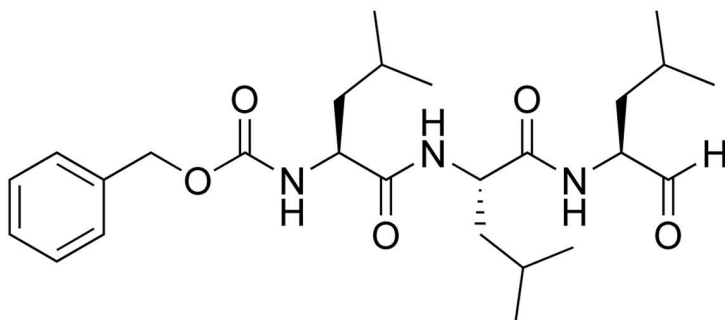

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com .
Solubility:	· DMSO ≤ 60 mM · Absolute ethanol ≤ 50 mM For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 210 μ L of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

CANCER RESEARCH

- Blocks apoptosis triggered by DNA damage in HeLa cells (Zhang et al.).
- Inhibits NF- κ B activation, sensitizing a variety of carcinoma cell lines to apoptosis (Arlt et al.).
- Cytotoxic effects on a variety of human cancer cell lines (Banerjee & Liefshitz).
- Inhibits growth of mouse melanoma (B16) and human ocular melanoma (IPC227F) cell lines (Vivier et al.).

References

- Arlt A et al. (2001) Inhibition of NF- κ B sensitizes human pancreatic carcinoma cells to apoptosis induced by etoposide (VP16) or doxorubicin. *Oncogene* 20(7): 859–68.
- Banerjee D & Liefshitz A. Potential of the proteasomal inhibitor MG-132 as an anticancer agent, alone and in combination. *Anticancer Res* 21(6A): 3941–7.
- Kisselev AF & Goldberg AL. (2001) Proteasome inhibitors: from research tools to drug candidates. *Chem Biol* 8(8): 739–58.
- Ortiz-Lazareno PC et al. (2008) MG132 proteasome inhibitor modulates proinflammatory cytokines production and expression of their receptors in U937 cells: involvement of nuclear factor-kappaB and activator protein-1. *Immunology* 124(4): 534–41.
- Palombella VJ et al. (1994) The ubiquitin-proteasome pathway is required for processing the NF-kappa B1 precursor protein and the activation of NF-kappa B. *Cell* 78(5): 773–85.
- Tsubuki S et al. (1996) Differential inhibition of calpain and proteasome activities by peptidyl aldehydes of di-leucine and tri-leucine. *J Biochem* 119(3): 572–6.
- Vivier M et al. (2008) Synthesis, radiosynthesis, and biological evaluation of new proteasome inhibitors in a tumor targeting approach. *J Med Chem* 51(4): 1043–7.
- Zhang L et al. (2011) MG132 inhibition of proteasome blocks apoptosis induced by severe DNA damage. *Cell Cycle* 10(20): 3515–8.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.