

Small Molecules

MK0457

Aurora kinase pathway inhibitor;
Inhibits aurora kinase A

Catalog # 73282
73284

50 mg
250 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM


FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

MK0457 is a pan-Aurora kinase inhibitor with a preference for Aurora A ($K_i = 0.6$ nM) over Aurora B ($K_i = 18$ nM) or Aurora C ($K_i = 4.6$ nM; Pollard & Mortimore; Harrington et al.). It is selective in a protein kinase screen (190 kinases tested) with the closest hit being Fms-related tyrosine kinase-3 (FLT-3; $K_i = 30$ nM; Pollard & Mortimore; Harrington et al.).

Molecular Name: MK0457
Alternative Names: Tozasertib, VX 680
CAS Number: 639089-54-6
Chemical Formula: $C_{23}H_{28}N_8OS$
Molecular Weight: 464.6 g/mol
Purity: $\geq 98\%$
Chemical Name: N-[4-[4-(4-methylpiperazin-1-yl)-6-[(5-methyl-1H-pyrazol-3-yl)amino]pyrimidin-2-yl]sulfanylphenyl]cyclopropanecarboxamide

Structure:

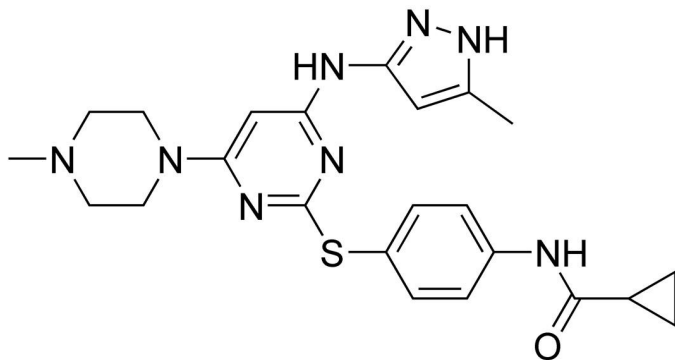

Properties

Physical Appearance: A crystalline solid
Storage: Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com.
Solubility: $\cdot$ DMSO ≤ 30 mM
For example, to prepare a 10 mM stock solution in DMSO, resuspend 50 mg in 10.8 mL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C . Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

CANCER RESEARCH

- Inhibits proliferation of a wide variety of tumor cell lines in vitro (Harrington et al.; Li et al.; Arlot-Bonnemains et al.).
- Blocks the growth of tumors in xenograft models of cancer, inhibits histone H3 phosphorylation, and increases apoptosis (Harrington et al.; Li et al.; Lin et al.).
- Disrupts bipolar spindle formation during mitosis, arresting cell cycle progression at the G2/M phase (Li et al.).

References

- Arlot-Bonnemains Y et al. (2008) Effects of the Aurora kinase inhibitor VX-680 on anaplastic thyroid cancer-derived cell lines. *Endocr Relat Cancer* 15(2): 559–68.
- Harrington EA et al. (2004) VX-680, a potent and selective small-molecule inhibitor of the Aurora kinases, suppresses tumor growth in vivo. *Nat Med* 10(3): 262–7.
- Li Y et al. (2010) VX680/MK-0457, a potent and selective Aurora kinase inhibitor, targets both tumor and endothelial cells in clear cell renal cell carcinoma. *Am J Transl Res* 2(3): 296–308.
- Lin YG et al. (2008) Targeting aurora kinase with MK-0457 inhibits ovarian cancer growth. *Clin Cancer Res* 14(17): 5437–46.
- Pollard JR & Mortimore M. (2009) Discovery and development of aurora kinase inhibitors as anticancer agents. *J Med Chem* 52(9): 2629–51.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.