

Small Molecules

Paclitaxel

Inhibitor of microtubule formation

Catalog # 73312
73314

25 mg
100 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM


FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Paclitaxel is a diterpene alkaloid originally isolated from the bark of the Pacific Yew tree (*Taxus brevifolia*). It binds to and stabilizes microtubules, preventing their reorganization during cell division, which leads to cell cycle arrest. Paclitaxel has antitumorigenic properties and has been used as a chemotherapeutic compound (Rowinsky et al.). Many pathways have been implicated in Paclitaxel-induced apoptosis, including c-Jun N-terminal kinase/stress-activated protein kinase (JNK/SAPK), p38 mitogen-activated protein kinase (MAPK), and protein kinase A (PKA; Wang et al.; Reshkin et al.).

Molecular Name: Paclitaxel
Alternative Names: NSC 125973
CAS Number: 33069-62-4
Chemical Formula: C₄₇H₅₁NO₁₄
Molecular Weight: 853.9 g/mol
Purity: ≥ 98%
Chemical Name: (2aR,4S,4aS,6R,9S,11S,12S,12aR,12bS)-1,2a,3,4,4a,6,9,10,11,12,12a,12b-Dodecahydro-4,6,9,11,12,12b-hexahydroxy-4a,8,13,13-tetramethyl-7,11-methano-5H-cyclodeca(3,4)benz(1,2-b)oxet-5-one 6,12b-diacetate, 12-benzoate, 9-ester with (2R,3S)-N-benzoyl-3-phenylisoserine

Structure:

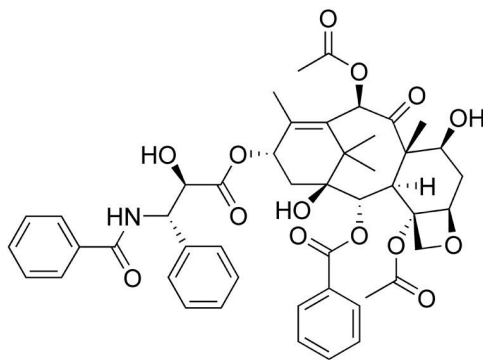

Properties

Physical Appearance: A crystalline solid
Storage: Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com.
Solubility:

- DMSO ≤ 5 mM
- Absolute ethanol ≤ 1.5 mM

For example, to prepare a 1 mM stock solution in DMSO, resuspend 10 mg in 11.7 mL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

DIFFERENTIATION

- Inhibits initiation and outgrowth of neurites in vitro, through microtubule polymerization (Letourneau & Ressler).

CANCER RESEARCH

- Inhibits tumor cell growth in a variety of cancer cell lines including cervical (HeLa), lung (A549), breast (MCF-7), colon (HT-29), ovarian (OVG-1), and pancreatic (PC-Sh) carcinomas (Liebmann et al.).
- Induces abnormal multipolar spindle formation, inducing cell cycle arrest at prophase and G1 in various human cell cancer lines (Woods et al.).
- Initiates apoptosis of cancer cells through multiple mechanisms involving: p53-dependent and -independent pathways, B-cell CLL/lymphoma 2 (BCL-2) family members, cyclin-dependent kinases, p38 MAPK, PKA, and JNK/SAPK (Wang et al.; Reshkin et al.).
- Induces cyclin inhibitor p21 in MCF7 and PC3M human cancer cell lines by a mechanism dependent on the activation of RAF-1 (Blagosklonny et al.).

References

- Blagosklonny M V et al. (1995) Taxol induction of p21WAF1 and p53 requires c-raf-1. *Cancer Res* 55(20): 4623–6.
- Letourneau PC & Ressler AH. (1984) Inhibition of neurite initiation and growth by taxol. *J Cell Biol* 98(4): 1355–62.
- Liebmann JE et al. (1993) Cytotoxic studies of paclitaxel (Taxol) in human tumour cell lines. *Br J Cancer* 68(6): 1104–9.
- Reshkin SJ et al. (2003) Paclitaxel induces apoptosis via protein kinase A- and p38 mitogen-activated protein-dependent inhibition of the Na⁺/H⁺ exchanger (NHE) NHE isoform 1 in human breast cancer cells. *Clin Cancer Res* 9(6): 2366–73.
- Rowinsky EK et al. (1990) Taxol: a novel investigational antimicrotubule agent. *J Natl Cancer Inst* 82(15): 1247–59.
- Wang TH et al. (2000) Paclitaxel-induced cell death: where the cell cycle and apoptosis come together. *Cancer* 88(11): 2619–28.
- Woods CM et al. (1995) Taxol-induced mitotic block triggers rapid onset of a p53-independent apoptotic pathway. *Mol Med* 1(5): 506–26.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.