

Small Molecules

Puromycin

Antibiotic; Protein synthesis inhibitor

Catalog # 73342
73344

50 mg
500 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM


FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Puromycin is an aminonucleoside antibiotic, derived from *Streptomyces alboniger*, which acts as a protein synthesis inhibitor. It binds to the target ribosome site A, where it is transferred to the growing polypeptide chain causing premature chain termination (Lührmann et al.; Rodriguez-Fonseca et al.; Azzam & Algranati). A biologically inactive form is generated when puromycin is N-acetylated by puromycin-N-acetyltransferase, allowing this gene to be used as a selective resistance marker (Vara et al.).

Molecular Name:	Puromycin (Dihydrochloride)
Alternative Names:	CL 13900, CL 16536, NSC 3055, PDH
CAS Number:	58-58-2
Chemical Formula:	$C_{22}H_{29}N_7O_5 \cdot 2HCl$
Molecular Weight:	544.4 g/mol
Purity:	≥ 98%
Chemical Name:	(2S)-2-amino-N-[(2S,4R,5R)-5-[6-(dimethylamino)purin-9-yl]-4-hydroxy-2-(hydroxymethyl)oxolan-3-yl]-3-(4-methoxyphenyl)propanamide;dihydrochloride

Structure:

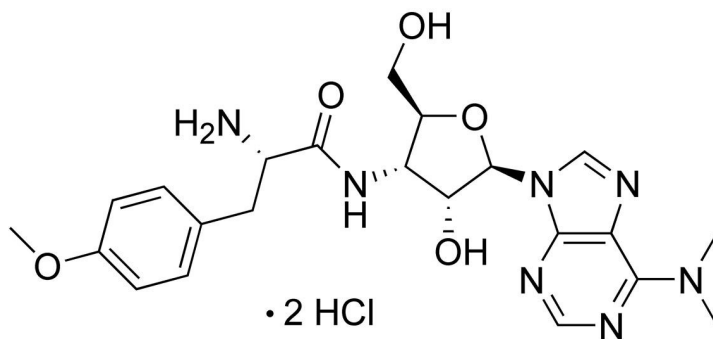

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com .
Solubility:	<ul style="list-style-type: none">· PBS (pH 7.2) ≤ 15 mM· DMSO ≤ 20 mM· Absolute ethanol ≤ 1.5 mM For example, to prepare a 10 mM stock solution in PBS, resuspend 50 mg in 9.18 mL of PBS (pH 7.2).

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

CELL LINE DEVELOPMENT

- Selects for cells expressing puromycin-N-acetyltransferase resistance gene as a research tool (de la Luna & Ortín; Iwaki et al.).
- Useful in CRISPR/Cas9 mammalian gene editing by selecting for successful Cas9-induced knock-in with puromycin resistance gene (Park et al.).

CANCER RESEARCH

- Possesses anti-tumor activity when tested against numerous cell lines (Foley & Eagle).

References

- Azzam ME & Algranati ID. (1973) Mechanism of puromycin action: fate of ribosomes after release of nascent protein chains from polysomes. *Proc Natl Acad Sci U S A* 70(12): 3866–9.
- Foley GE & Eagle H. (1958) The cytotoxicity of anti-tumor agents for normal human and animal cells in first tissue culture passage. *Cancer Res* 18(9): 1011–6.
- Iwaki T et al. (2003) Rapid selection of *Drosophila* S2 cells with the puromycin resistance gene. *Biotechniques* 35(3): 482–4, 486.
- de la Luna S & Ortín J. (1992) Recombinant DNA Part G Elsevier.
- Lührmann R et al. (1981) Localization of the puromycin binding site on the large ribosomal subunit of *Escherichia coli* by immunoelectron microscopy. *Proc Natl Acad Sci U S A* 78(12): 7276–80.
- Park A et al. (2014) CRISPR/Cas9 allows efficient and complete knock-in of a destabilization domain-tagged essential protein in a human cell line, allowing rapid knockdown of protein function. *PLoS One* 9(4): e95101.
- Rodriguez-Fonseca C et al. (2000) Puromycin-rRNA interaction sites at the peptidyl transferase center. *RNA* 6(5): 744–54.
- Vara et al. (1985) Biosynthesis of puromycin by *Streptomyces alboniger*: Characterization of puromycin N-acetyltransferase. *Biochemistry* 24(27): 8074–8081.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.