

Small Molecules

SU11274

HGF pathway inhibitor; Inhibits MET



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog # 73432
73434


1 mg
5 mg

Product Description

SU11274 is a selective, ATP-competitive inhibitor of c-MET receptor tyrosine kinase ($IC_{50} = 20$ nM; Sattler et al.) It shows good selectivity towards c-MET versus other receptor tyrosine kinases with IC_{50} values of $1.3 \mu\text{M}$ and $9.7 \mu\text{M}$ for fetal liver kinase (FLK-1) and fibroblast growth factor receptor-1 (FGFR-1), respectively (Sattler et al.). It retains activity for certain c-MET mutants, for example H1112Y and M1268T (Berthou et al.).

Molecular Name: SU11274
Alternative Names: Met Kinase Inhibitor
CAS Number: 658084-23-2
Chemical Formula: $C_{28}H_{30}ClN_5O_4S$
Molecular Weight: 568.1 g/mol
Purity: $\geq 98\%$
Chemical Name: (3Z)-N-(3-chlorophenyl)-3-[[[3,5-dimethyl-4-[(4-methyl-1-piperazinyl)carbonyl]-1H-pyrrol-2-yl]methylene]-2,3-dihydro-N-methyl-2-oxo-1H-indole-5-sulfonamide

Structure:

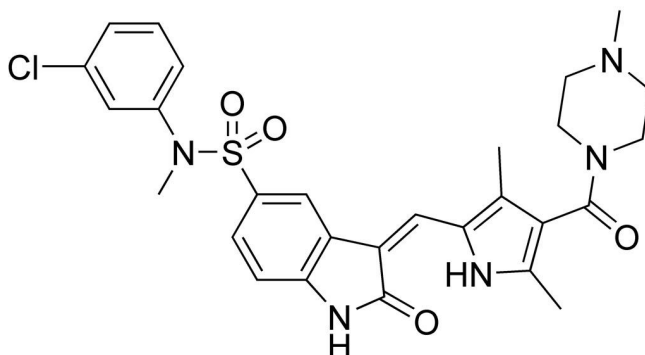

Properties

Physical Appearance: A crystalline solid
Storage: Product stable at -20°C as supplied. Protect from prolonged exposure to light. For product expiry date, please contact techsupport@stemcell.com.
Solubility: $\cdot$ DMSO ≤ 15 mM
 $\cdot$ Absolute ethanol ≤ 0.4 mM
For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 176 μL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C . Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

CANCER RESEARCH

- Induces apoptosis and cell cycle arrest in TPR-MET-transformed mouse pre-B (Ba/F3) cells (Sattler et al.).
- Inhibits cell viability and motility of human head and neck squamous cell carcinoma cell lines in vitro (Seiwert et al.).
- Inhibits cell viability of c-MET-expressing human non-small cell lung carcinoma cell lines in vitro (Ma et al.).

References

- Berthou S et al. (2004) The Met kinase inhibitor SU11274 exhibits a selective inhibition pattern toward different receptor mutated variants. *Oncogene* 23(31): 5387–93.
- Ma et al. (2005) Functional Expression and Mutations of c-Met and its therapeutic inhibition with SU11274 and small interfering RNA in non-small cell lung cancer. *Cancer Res* 65(4): 1479-88.
- Sattler M et al. (2003) A novel small molecule met inhibitor induces apoptosis in cells transformed by the oncogenic TPR-MET tyrosine kinase. *Cancer Res* 63(17): 5462–9.
- Seiwert TY et al. (2009) The MET receptor tyrosine kinase is a potential novel therapeutic target for head and neck squamous cell carcinoma. *Cancer Res* 69(7): 3021–31.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.