

Cytokines

Human Recombinant HGF



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Hepatocyte growth factor

Catalog #	78019.1	10 µg
	78019	50 µg
	78019.2	1000 µg

Product Description

Hepatocyte growth factor/scatter factor (HGF/SF) is a pleiotropic growth factor that promotes proliferation, motility, survival, and differentiation. HGF binds and promotes the dimerization and activation of the receptor tyrosine kinase c-MET, and stimulates PI3K/AKT, FAK, JNK, and ERK 1/2 signaling (Organ & Tsao). HGF stimulates migration of cells during embryogenesis, induces cell motility and scattering of epithelial cells, and regulates epithelial-mesenchymal transition. It is a chemoattractant for motor neurons, and regulates the development of sensory, sympathetic, parasympathetic, and cortical neurons as well as the proliferation of oligodendrocytes and glial development. Additionally, HGF is cytoprotective and regulates liver growth and regeneration (Nakamura et al.).

Product Information

Alternative Names:	DFNB39, F-TCF, Fibroblast-derived tumor cytotoxic factor, Hepapoiectin A, HPTA, Lung fibroblast-derived mitogen, Scatter factor, SF
Accession Number:	P14210
Amino Acid Sequence:	QRKRRNTIHE FKKSATTLI KIDPALKIKT KKVNTADQCA NRCTRNGKLP FTCKAFVFDK ARKQCLWFPF NSMSSGVKKE FGHEFDLYEN KDYIRNCIIG KGRSYKGTVS ITKSGIKCQP WSSMIPHEHS FLPSSYRGKD LQENYCRNPR GEEGPPWCFT SNPEVRYEVC DIPQCSEVEC MTCNGESYRG LMDHTESGKI QQRWDHQTPH RHKFLPERYP DKGFDNYCR NPDGQPRPWC YTLDPHTRWE YCAIKTCADN TMNDTDVPLE TTECIQGQGE GYRGTVNTIW NGIPCQRWDS QYPHEHDMTP ENFKCKDLRE NYCRNPDGSE SPWCFTTDPN IRVGYCSQIP NCDMSHGQDC YRNGGKNYMG NLSQTRSGLT CSMWDKNMED LHRHIFWEPD ASKLNENYCR NPDDDAHGPW CYTGNPLIPW DYCPISRCEG DTTPTIVNLD HPVISCAKTK QLRVVNGIPT RTNIGWMVSL RYRNKHICGG SLIKESWWLT ARQCFSRDL KDYEAWLGIH DVHGRGDEKC KQVLNVSQVL YGPEGSDLVL MKLARPAVLD DFVSTIDLPN YGCTIPEKTS CSVYGWGYTG LINYDGLLRV AHLYIMGNEK CSQHHRGKVT LNESEICAGA EKIGSGPCEG DYGGPLVCEQ HKMRMVLGVI VPGRGCAIPN RPGIFVRVAY YAKWIHKIIL TYKVPQS
Predicted Molecular Mass:	88 - 90 kDa
Species:	Human
Cross Reactivity:	Mouse, Rat
Formulation:	Lyophilized after dialysis against phosphate-buffered saline.
Source:	CHO


Specifications

Activity:	The specific activity is $\geq 1 \times 10^5$ units/mg ($EC_{50} \leq 10$ ng/mL) as determined by a cell proliferation assay using 4MBr5 cells.
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic Limulus amoebocyte lysate (LAL) analysis and is ≤ 0.2 EU/µg protein.

Preparation and Storage

Storage:	Store at -80°C.
Stability:	Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Resuspend the product in sterile water or phosphate-buffered saline to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. Store at 2 - 8°C for up to 1 week or at -20°C to -80°C for up to 2 months. Avoid repeated freeze-thaw cycles.

Data

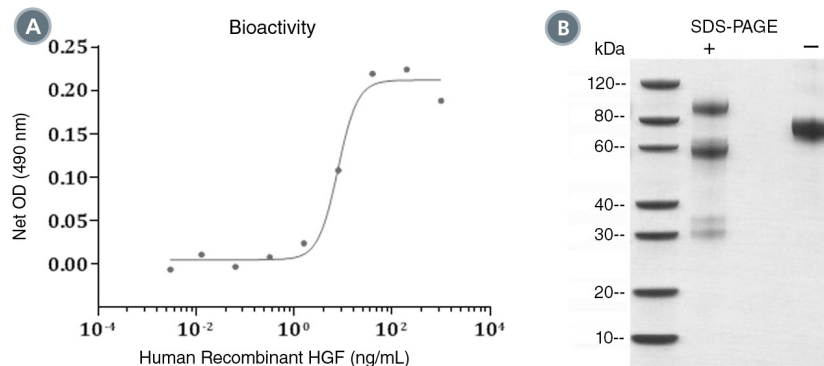

(A) The biological activity of Human Recombinant HGF was tested by its ability to promote the proliferation of 4MBr5 cells. Cell proliferation was measured using a fluorometric assay method. The EC50 is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC50 in the above example is 7.2 ng/mL.

(B) 5 µg of Human Recombinant HGF was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Human Recombinant HGF has a predicted molecular mass of 88 - 90 kDa (single chain), 59 - 61 kDa (alpha chain), or 30 - 34 kDa (beta chain).

Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Nakamura T et al. (2011) Hepatocyte growth factor twenty years on: Much more than a growth factor. *J Gastroenterol Hepatol* 26 Suppl 1: 188–202.
- Organ SL & Tsao M-S. (2011) An overview of the c-MET signaling pathway. *Ther Adv Med Oncol* 3(1 Suppl): S7–S19.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.