

## Cytokines

### Mouse Recombinant TNF-alpha

Tumor necrosis factor alpha



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog # 78069.1  
78069

20 µg  
100 µg

## Product Description

Tumor necrosis factor- $\alpha$  (TNF- $\alpha$ ) is a pro-inflammatory cytokine that activates NF- $\kappa$ B, MAPK, and PI3K/AKT pathways. Activated T cells and macrophages are the primary producers of TNF- $\alpha$  in response to inflammation and infectious conditions. Many other cell types have been shown to produce TNF- $\alpha$ , among them B cells, NK cells, mast cells, neutrophils, dendritic cells, microglia, endothelial cells, smooth muscle cells, cardiomyocytes, and fibroblasts. TNF- $\alpha$  has cytotoxic effects on cancer cells in vitro by stimulating anti-tumor immunosuppressive responses. TNF- $\alpha$  stimulates expression of E- and P-selectins, thus facilitating adhesion of neutrophils, monocytes, and memory T cells to activated platelets and endothelial cells (Zelová & Hosek). Other effects of TNF- $\alpha$  include vasodilatation and edema formation.

## Product Information

**Alternative Names:** Cachectin, Cachexin, Cytotoxin, DIF, Necrosin, TNF, TNF- $\alpha$ , TNFSF2, Tumor necrosis factor- $\alpha$   
**Accession Number:** P06804  
**Amino Acid Sequence:** MLRSSSQNSS DKPVAHVVAN HQVEEQLEWL SQRANALLAN GMDLKDNQLV VPADGLYLVY SQVLFKGQGC PDYVLLTHTV SRFAISYQEK VNLLSAVKSP CPKDTPEGAE LKPWYEPIYL GGVFQLEKGD QLSAEVNLPK YLDAESGQV YFGVIAL  
**Predicted Molecular Mass:** 17.4 kDa  
**Species:** Mouse  
**Cross Reactivity:** Human, Rat, Monkey  
**Formulation:** Lyophilized from a sterile-filtered aqueous solution containing sodium phosphate and sodium chloride, pH 7.5.  
**Source:** E. coli

## Specifications


**Activity:** The specific activity is  $\geq 1 \times 10^7$  units/mg ( $EC_{50} \leq 0.1$  ng/mL) as determined by the cytolysis of mouse L929 cells growing in the presence of actinomycin D.  
**Purity:**  $\geq 95\%$   
**Endotoxin Level:** Measured by kinetic Limulus amebocyte lysate (LAL) analysis and is  $\leq 1$  EU/µg protein.

## Preparation and Storage

**Storage:** Store at -20°C to -80°C.  
**Stability:** Stable as supplied for 12 months from date of receipt.  
**Preparation:** Centrifuge vial before opening. Reconstitute the product in sterile water to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex.

OPTIONAL: After reconstitution, if product will not be used immediately, dilute with concentrated bovine serum albumin (BSA) to a final BSA concentration of 0.1%. The effect of storage of stock solution on product performance should be tested for each application. As a general guide, do not store at 2 - 8°C for more than 1 month or at -80°C for more than 3 months. Avoid repeated freeze-thaw cycles.

## Data

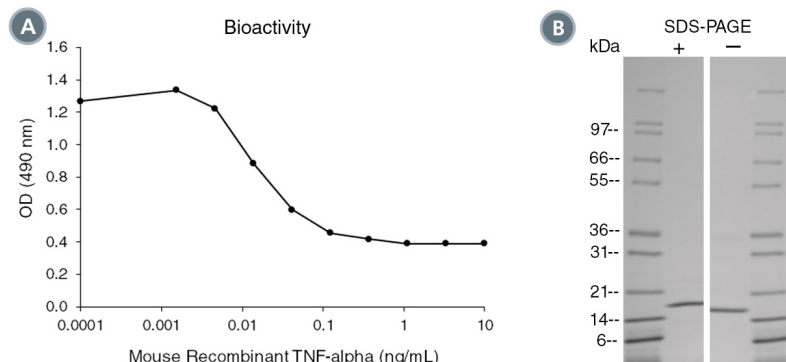

(A) The biological activity of Mouse Recombinant TNF-alpha was tested by its ability to induce cytolysis of mouse L929 cells growing in the presence of actinomycin D. Cell viability was measured after 19 hours of culture using a fluorometric assay method. The EC50 is defined as the effective concentration of the growth factor at which cell viability is at 50% of maximum. The EC50 in the above example is 0.013 - 0.02 ng/mL.

(B) 1 µg of Mouse Recombinant TNF-alpha was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Mouse Recombinant TNF-alpha has a predicted molecular mass of 17.4 kDa.

## Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, visit [www.stemcell.com/cytokines](http://www.stemcell.com/cytokines) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

## References

Zelová H & Hošek J. (2013) TNF- $\alpha$  signalling and inflammation: interactions between old acquaintances. *Inflamm Res* 62(7): 641–51.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2018 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.