

## Anti-Human CD4 Antibody, Clone SK3, PE



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Antibodies

Mouse monoclonal IgG1 antibody  
against human, rhesus, cynomolgus  
CD4, PE-conjugated

Catalog #60122PE  
#60122PE.1

100 Tests 5 µL/test  
25 Tests 5 µL/test

## Product Description

The SK3 antibody reacts with CD4, an ~59 kDa single-chain type 1 transmembrane glycoprotein and member of the immunoglobulin (Ig) superfamily; CD4 contains four extracellular Ig-like domains (D1 - D4). The epitope for SK3 has been localized to the D3 domain of the protein, which has a structure resembling an Ig variable domain. CD4 is expressed at relatively high levels by most thymocytes and a subpopulation of T cells (T-helper cells), and at lower levels by peripheral blood monocytes and macrophages. CD4 binds to a nonpolymorphic region of MHC II and acts as a co-receptor to the T cell receptor (TCR) in MHC II-restricted antigen recognition by enhancing the affinity of the association between the TCR and MHC II-antigen complex. CD4 also functions to amplify signals from the TCR to the cytoplasm through the interaction of its intracellular domain with cytoplasmic tyrosine kinases such as Lck. Moreover, CD4 is a receptor for the human immunodeficiency virus (HIV). Binding of the SK3 antibody blocks HIV binding and mixed lymphocyte reaction (MLR).

|                      |                                                                                                     |
|----------------------|-----------------------------------------------------------------------------------------------------|
| Target Antigen Name: | CD4                                                                                                 |
| Alternative Names:   | Leu3, T4                                                                                            |
| Gene ID:             | 920                                                                                                 |
| Species Reactivity:  | Human, Rhesus, Cynomolgus, African green, Olive baboon, Pigtailed macaque, Sooty mangabey, Squirrel |
| Host Species:        | Mouse (BALB/c)                                                                                      |
| Clonality:           | Monoclonal                                                                                          |
| Clone:               | SK3                                                                                                 |
| Isotype:             | IgG1, kappa                                                                                         |
| Immunogen:           | Human peripheral blood T lymphocytes                                                                |
| Conjugate:           | PE                                                                                                  |

## Applications


|           |    |
|-----------|----|
| Verified: | FC |
| Reported: | FC |

Abbreviations: CellSep: Cell separation; ChIP: Chromatin immunoprecipitation; FA: Functional assay; FC: Flow cytometry; ICC: Immunocytochemistry; IF: Immunofluorescence microscopy; IHC: Immunohistochemistry; IP: Immunoprecipitation; RIA: Radioimmunoassay; WB: Western blotting

## Properties

|                        |                                                                                                                                                                                                                                          |
|------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Formulation:           | Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) bovine serum albumin                                                                                                                                   |
| Purification:          | The antibody was purified by affinity chromatography and conjugated with PE under optimal conditions. The solution is free of unconjugated PE and unconjugated antibody.                                                                 |
| Stability and Storage: | Product stable at 2 - 8°C when stored undiluted. Do not freeze. Protect product from prolonged exposure to light. For product expiry date, please contact <a href="mailto:techsupport@stemcell.com">techsupport@stemcell.com</a> .       |
| Directions for Use:    | For flow cytometry the suggested use of this antibody is 5 µL per 1 x 10 <sup>6</sup> cells in 100 µL volume or per 100 µL of whole blood. It is recommended that the antibody be titrated for optimal performance for each application. |

## Data

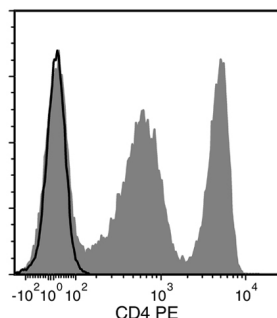

Flow cytometry analysis of human peripheral blood mononuclear cells (PBMCs) labeled with Anti-Human CD4 Antibody, Clone SK3, PE (filled histogram) or Mouse IgG1, kappa Isotype Control Antibody, Clone MOPC-21, PE (Catalog #60070PE; solid line histogram).

## Related Products

For a complete list of antibodies, including other conjugates, sizes and clones, as well as related products available from STEMCELL Technologies, please visit our website at [www.stemcell.com/antibodies](http://www.stemcell.com/antibodies) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

## References

1. Wang L et al. (2012) Human CD4+ lymphocytes for antigen quantification: characterization using conventional flow cytometry and mass cytometry. *Cytometry A* 81(7): 567–75. (CyTOF, FC)
2. Green AM et al. (2010) CD4+ regulatory T cells in a cynomolgus macaque model of *Mycobacterium tuberculosis* infection. *J Infect Dis* 202(4): 533–41. (FC)
3. Chen X et al. (2009) Engagement of Toll-like receptor 2 on CD4+ T cells facilitates local immune responses in patients with tuberculous pleurisy. *J Infect Dis* 200(3): 399–408. (FC)
4. Basch RS et al. (1990) Expression of CD4 by human megakaryocytes. *Proc Natl Acad Sci USA* 87(20): 8085–9. (FC, IF)
5. McDougal JS et al. (1986) Binding of the human retrovirus HTLV-III/LAV/ARV/HIV to the CD4 (T4) molecule: Conformation dependence, epitope mapping, antibody inhibition, and potential for idiotypic mimicry. *J Immunol* 137(9): 2937–44. (IF, WB)
6. Sattentau QJ et al. (1986) Epitopes of the CD4 antigen and HIV infection. *Science* 234(4780): 1120–3. (FA)
7. Engleman EG et al. (1981) Antibodies to membrane structures that distinguish suppressor/cytotoxic and helper T lymphocyte subpopulations block the mixed leukocyte reaction in man. *J Exp Med* 154(1): 193–8. (FA, FC, Isolation of T Cell Subsets using a panning technique)
8. Ledbetter JA et al. (1981) Evolutionary conservation of surface molecules that distinguish T lymphocyte helper/inducer and cytotoxic/suppressor subpopulations in mouse and man. *J Exp Med* 153(2): 310–23. (FC, IP)

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.