

Cytokines

Human Recombinant BDNF

Brain-derived neurotrophic factor



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #	78005	10 µg
	78005.1	100 µg
	78005.2	1000 µg

Product Description

Brain-derived neurotrophic factor (BDNF), like nerve growth factor (NGF), neurotrophin-3 (NT-3), and neurotrophin-4 (NT-4), is a member of the NGF family of neurotrophins, which is required for the differentiation and survival of specific neuronal subpopulations in both the central and the peripheral nervous systems (Minichiello & Klein; Minichiello et al.). BDNF binds with high affinity to the TRKB kinase receptor, and activates AKT and ERK pathways (Mattson et al.). It is expressed in hippocampus, cortex, and synapses of the basal forebrain. BDNF acts as a survival factor for human embryonic stem cells when plated on either feeders or Corning® Matrigel® (Pyle et al.). BDNF regulates synaptic transmission and plasticity at adult synapses in the central nervous system, contributes to adaptive neuronal responses including long-term potentiation, long-term depression, certain forms of short-term synaptic plasticity, as well as homeostatic regulation of neuronal excitability (Reichardt). It also has a role in neurogenesis by promoting survival and growth of dorsal root ganglion cells, hippocampal, and cortical neurons (Binder & Scharfman). BDNF, together with Human Recombinant GDNF (glial cell line-derived neurotrophic factor; Catalog #78058), BrainPhys™ Neuronal Medium (Catalog #05790), and other supplements, can be used to differentiate human pluripotent stem cell (hPSC)-derived neural progenitor cells into neurons (Bardy et al.).

Product Information

Alternative Names:	Abrineurin, ANON2, BULN2, Neurotrophin, MGC34632
Accession Number:	P23560
Amino Acid Sequence:	MHSDPARRGE LSVCDISEW VTAADKKTAV DMSGGTVTVL EKVPVSKGQL KQYFYETKCN PMGYTKEGCR GIDKRHWNSQ CRTTQSYVRA LTMDSKKRIG WRFIRIDTSC VCTLTIKRGR
Predicted Molecular Mass:	13.6 kDa monomer; 27.3 kDa dimer
Species:	Human
Cross Reactivity:	Mouse, Rat
Formulation:	Lyophilized from a sterile filtered solution containing 0.1% trifluoroacetic acid.
Source:	E. coli


Specifications

Activity:	The specific activity is $\geq 5 \times 10^2$ units/mg ($EC_{50} \leq 2 \mu\text{g/mL}$) as determined by a cell proliferation assay of C6 cells.
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic limulus amoebocyte lysate (LAL) analysis and is ≤ 1 EU/ μg protein.

Preparation and Storage

Storage:	Store at -20°C to -80°C.
Stability:	Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Resuspend the product in sterile water containing 0.1% bovine serum albumin (BSA) to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. Store at 2 - 8°C for up to 1 month or at -20°C to -80°C for up to 3 months. Avoid repeated freeze-thaw cycles. NOTE: If reconstituted product will be used immediately BSA is not required.

Data

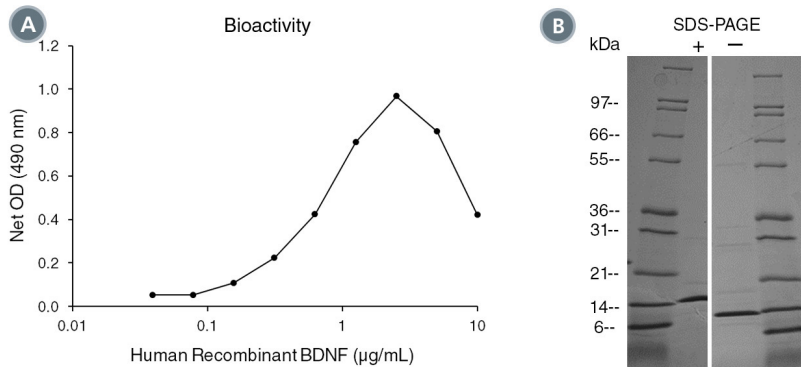

(A) The biological activity of Human Recombinant BDNF was tested by its ability to promote the proliferation of C6 cells. Cell proliferation was measured after 7 days of culture using a fluorometric assay method. The EC50 is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC50 in the example above is 0.65 µg/mL.

(B) 1 µg of Human Recombinant BDNF was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Human Recombinant BDNF is a homodimer of 13.6 kDa subunits with a predicted total molecular mass of 27.3 kDa.

Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Bardy C et al. (2015) Neuronal medium that supports basic synaptic functions and activity of human neurons in vitro. *Proc Natl Acad Sci U S A* 112(20): E2725–34.
- Binder DK & Scharfman HE. (2004) Brain-derived neurotrophic factor. *Growth Factors* 22(3): 123–31.
- Mattson MP et al. (2004) A neural signaling triumvirate that influences ageing and age-related disease: insulin/IGF-1, BDNF and serotonin. *Ageing Res Rev* 3(4): 445–64.
- Minichiello L & Klein R. (1996) TrkB and TrkC neurotrophin receptors cooperate in promoting survival of hippocampal and cerebellar granule neurons. *Genes Dev* 10(22): 2849–58.
- Minichiello L et al. (1995) Differential effects of combined trk receptor mutations on dorsal root ganglion and inner ear sensory neurons. *Development* 121(12): 4067–75.
- Pyle AD et al. (2006) Neurotrophins mediate human embryonic stem cell survival. *Nat Biotechnol* 24(3): 344–50.
- Reichardt LF. (2006) Neurotrophin-regulated signalling pathways. *Philos Trans R Soc Lond B Biol Sci* 361(1473): 1545–64.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2016 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.