

## Cytokines

### Human Recombinant IL-9



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #	78085	10 µg
	78085.1	50 µg
	78085.2	100 µg
	78085.3	500 µg
	78085.4	1000 µg

## Product Description

Interleukin 9 (IL-9) has pleiotropic functions in the immune system and signals through its specific IL-9 receptor (IL-9R) (Renauld et al.). IL-9/IL-9R signaling pathway mainly targets the downstream activation of JAK/STAT (janus kinase/signal transducer and activator of transcription) and subsequent phosphorylation cascades initiated by multiple kinases including IRS-PI3K-PKB (insulin receptor substrate, phosphatidylinositol 3-kinases, protein kinase-B) and ERK (Knoops & Renauld; Fontaine et al.). IL-9 has been shown to have a role in Th1/Th17-mediated inflammation and in regulatory T cell responses (Singh et al.; Goswami & Kaplan). IL-9/IL-9R signaling pathway represents a novel endogenous anti-apoptotic mechanism for cortical neurons (Fontaine et al.). IL-9 secreting T cells, termed Th9 cells, contribute to both effective immunity and immunopathological disease, and have been shown to have a role in the treatment of allergic and autoimmune disease (Kaplan et al.).

## Product Information

Alternative Names:	Interleukin-9, Mast cell growth-enhancing activity, Mast cell growth factor, MCGF, MEA, p40 cytokine, TCGF-3, T-cell growth factor 3
Accession Number:	P15248
Amino Acid Sequence:	QGCPFLAGIL DINFLINKMQ EDPASKCHCS ANVTSCCLCLG IPSDNCTRPC FSERLSQMTN TTMQTRYPLI FSRVKKSVEV LKNNKCPYFS CEQPCNQTTA GNALFLKSL LEIFQKEKMR GMRGKI
Predicted Molecular Mass:	25 - 40 kDa
Species:	Human
Cross Reactivity:	Mouse
Formulation:	Lyophilized after dialysis against phosphate-buffered saline.
Source:	CHO


## Specifications

Activity:	The specific activity is $\geq 1 \times 10^6$ units/mg ( $EC_{50} \leq 1$ ng/mL) as determined by a cell proliferation assay using MO7e cells.
Purity:	$\geq 95$ %
Endotoxin Level:	Measured by kinetic Limulus amebocyte lysate (LAL) analysis and is $\leq 0.2$ EU/µg protein.

## Preparation and Storage

Storage:	Store at -80°C.
Stability:	Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Resuspend the product in sterile water or phosphate-buffered saline to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. Store at 2 - 8°C for up to 1 week or at -20°C to -80°C for up to 2 months. Avoid repeated freeze-thaw cycles.

## Data

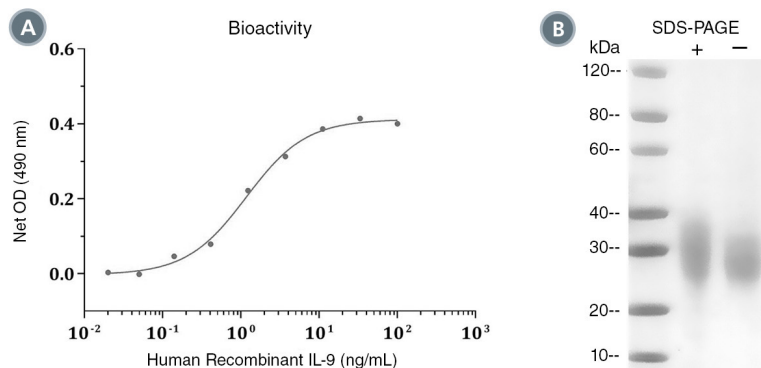

(A) The biological activity of Human Recombinant IL-9 was tested by its ability to promote the proliferation of MO7e cells. Cell proliferation was measured using a fluorometric assay method. The EC50 is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC50 in the above example is less than 1 ng/mL.

(B) 2 µg of Human Recombinant IL-9 was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining. Human Recombinant IL-9 has a predicted molecular mass of 25 - 40 kDa.

## Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, please visit our website at [www.stemcell.com/cytokines](http://www.stemcell.com/cytokines) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

## References

- Fontaine RH et al. (2008) IL-9/IL-9 receptor signaling selectively protects cortical neurons against developmental apoptosis. *Cell Death Differ* 15(10): 1542–52.
- Goswami R & Kaplan MH. (2011) A brief history of IL-9. *J Immunol* 186(6): 3283–8.
- Kaplan MH et al. (2015) The development and in vivo function of T helper 9 cells. *Nat Rev Immunol* 15(5): 295–307.
- Knoops L & Renauld J-C. (2004) IL-9 and its receptor: from signal transduction to tumorigenesis. *Growth Factors* 22(4): 207–15.
- Renauld JC et al. (1992) Expression cloning of the murine and human interleukin 9 receptor cDNAs. *Proc Natl Acad Sci U S A* 89(12): 5690–4.
- Singh TP et al. (2013) Involvement of IL-9 in Th17-associated inflammation and angiogenesis of psoriasis. *PLoS One* 8(1): e51752.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.