

Small Molecules

Dibutylryl-cAMP (Sodium Salt)

cAMP pathway activator; Activates cAMP-dependent protein kinases

Catalog #	73882	25 mg
	73884	100 mg
	73886	250 mg



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

Dibutylryl-cAMP is a cell-permeable cyclic AMP (cAMP) analog that activates cAMP-dependent protein kinases (Schwede et al.). This product is supplied as the sodium salt of the molecule.

Molecular Name:	Dibutylryl-cAMP (Sodium Salt)
Alternative Names:	Bucladesine; DC 2797
CAS Number:	16980-89-5
Chemical Formula:	$C_{18}H_{23}N_5O_8P \cdot Na$
Molecular Weight:	491.4 g/mol
Purity:	≥ 95%
Chemical Name:	N-(1-oxobutyl)-cyclic 3',5'-(hydrogen phosphate) 2'-butanoate-adenosine, monosodium salt
Structure:	

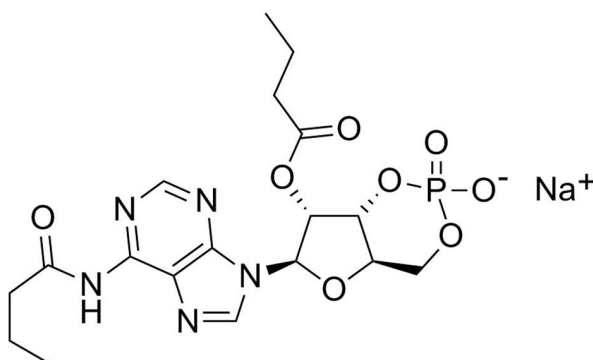

Properties

Physical Appearance:	A crystalline solid
Storage:	Product stable at -20°C as supplied. Protect product from prolonged exposure to light. For long-term storage store with a desiccant. For product expiry date, please contact techsupport@stemcell.com.
Solubility:	<ul style="list-style-type: none">· PBS ≤ 6.5 mM· DMSO ≤ 6.5 mM· Absolute ethanol ≤ 2 mM For example, to prepare a 5 mM stock solution in DMSO, resuspend 10 mg in 4.07 mL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

DIFFERENTIATION

- Suppresses experimental autoimmune encephalomyelitis development by reducing demyelination and mobilizing neural stem cells in the subventricular zone towards the demyelinated plaques (Khezri et al.).
- Induces intrinsic axon growth in peripheral and central nervous systems, and morphological differentiation of astrocytes (Knott et al.; Imamura et al.).
- Stimulates neurite outgrowth in PC12 cells (Maruoka et al.).

References

- Imamura M & Ozawa E. (1998) Differential expression of dystrophin isoforms and utrophin during dibutyryl-cAMP-induced morphological differentiation of rat brain astrocytes. *Proc Natl Acad Sci U S A* 95(11): 6139–44.
- Khezri S et al. (2013) Dibutyryl cyclic AMP inhibits the progression of experimental autoimmune encephalomyelitis and potentiates recruitment of endogenous neural stem cells. *J Mol Neurosci* 51(2): 298–306.
- Knott EP et al. (2014) Cyclic AMP signaling: a molecular determinant of peripheral nerve regeneration. *Biomed Res Int* 2014: Article ID 651625.
- Maruoka H et al. (2010) Dibutyryl-cAMP up-regulates nur77 expression via histone modification during neurite outgrowth in PC12 cells. *J Biochem* 148(1): 93–101.
- Schwede F et al. (2000) Cyclic nucleotide analogs as biochemical tools and prospective drugs. *Pharmacol Ther* 87(2–3): 199–226.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.