

Small Molecules

(±)-Flurbiprofen

PI3K/AKT pathway inhibitor; Inhibits COX-1 and COX-2

Catalog # 73922
73924

1 g
5 g



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Product Description

(±)-Flurbiprofen is a cell-permeable and non-steroidal anti-inflammatory agent. It is a potent non-selective inhibitor of cyclooxygenase (COX) with EC₅₀ = 0.1 and 0.4 μM for the inhibition of recombinant human COX-1 and COX-2, respectively (Gierse et al.).

Molecular Name: (±)-Flurbiprofen
Alternative Names: Ansaïd
CAS Number: 5104-49-4
Chemical Formula: C₁₅H₁₃FO₂
Molecular Weight: 244.3 g/mol
Purity: ≥ 99%
Chemical Name: 2-(3-fluoro-4-phenylphenyl)propanoic acid
Structure:

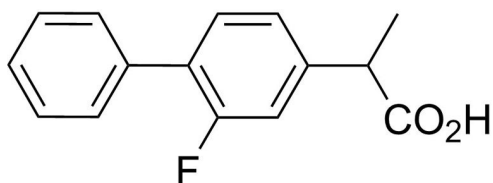

Properties

Physical Appearance: A crystalline solid
Storage: Product stable at room temperature (15 - 25°C) as supplied. Protect product from prolonged exposure to light. For long-term storage store with a desiccant. For product expiry date, please contact techsupport@stemcell.com.
Solubility: · DMSO ≤ 400 mM
· Absolute ethanol ≤ 400 mM
For example, to prepare a 10 mM stock solution in DMSO, resuspend 10 mg in 4.09 mL of DMSO.

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

Published Applications

MAINTENANCE AND SELF-RENEWAL

- Maintains self-renewal and pluripotency in hESC cells grown in conditions that induce differentiation (Bhanu et al.; Desbordes et al.).

CANCER RESEARCH

- Inhibits tumor growth in models of prostate cancer, familial adenomatous polyposis, and brain cancer (King, Jr & Khalili; Grubbs et al.; Wechter et al.; Zemskova et al.).

References

Bhanu N V. et al. (2016) Histone modification profiling reveals differential signatures associated with human embryonic stem cell self-renewal and differentiation. *Proteomics* 16(3): 448–58.

Desbordes SC et al. (2008) High-throughput screening assay for the identification of compounds regulating self-renewal and differentiation in human embryonic stem cells. *Cell Stem Cell* 2(6): 602–12.

Gierse JK et al. (1995) Expression and selective inhibition of the constitutive and inducible forms of human cyclo-oxygenase. *Biochem J*: 479–84.

Grubbs CJ et al. (2000) Celecoxib inhibits N-Butyl-N-(4-hydroxybutyl)-nitrosamine-induced urinary bladder cancers in male B6D2F1 mice and female Fischer-344 rats. *Cancer Res* 60(20): 5599-602.

King, Jr JG & Khalili K. (2001) Inhibition of human brain tumor cell growth by the anti-inflammatory drug, flurbiprofen. *Oncogene* 20(47): 6864–70.

Wechter WJ et al. (2000) E-7869 (R-Flurbiprofen) inhibits progression of prostate cancer in the TRAMP mouse. *Cancer Res* 60(8): 2203–8.

Zemskova M et al. (2006) Gene expression profiling in R-flurbiprofen-treated prostate cancer: R-Flurbiprofen regulates prostate stem cell antigen through activation of AKT kinase. *Biochem Pharmacol* 72(10): 1257–67.

Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at www.stemcell.com/smallmolecules or contact us at techsupport@stemcell.com.

This product is hazardous. Please refer to the Safety Data Sheet (SDS).

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.