

## Small Molecules

### Phorbol 12-myristate 13-acetate

PKC pathway activator; Activates protein kinase C

Catalog # 74042  
74044

1 mg  
5 mg



Scientists Helping Scientists™ | [WWW.STEMCELL.COM](http://WWW.STEMCELL.COM)

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713


[INFO@STEMCELL.COM](mailto:INFO@STEMCELL.COM) • [TECHSUPPORT@STEMCELL.COM](mailto:TECHSUPPORT@STEMCELL.COM)

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

## Product Description

Phorbol 12-myristate 13-acetate (PMA), a phorbol ester, is a reversible, highly potent activator of protein kinase C (PKC), including Group A and Group B isoforms (Blumberg). PMA activates PKC isoforms through the binding of the C1 domain (Hurley et al.). PMA can also activate certain mitogen-activated protein (MAP) kinase pathways through PKC (Adams & Parker).

|                    |                                                |
|--------------------|------------------------------------------------|
| Molecular Name:    | Phorbol 12-myristate 13-acetate                |
| Alternative Names: | 12-O-Tetradecanoylphorbol-13-acetate; PMA; TPA |
| CAS Number:        | 16561-29-8                                     |
| Chemical Formula:  | C <sub>36</sub> H <sub>56</sub> O <sub>8</sub> |
| Molecular Weight:  | 616.8 g/mol                                    |
| Purity:            | ≥ 98%                                          |
| Chemical Name:     | 12-O-tetradecanoylphorbol-13-acetate           |
| Structure:         |                                                |

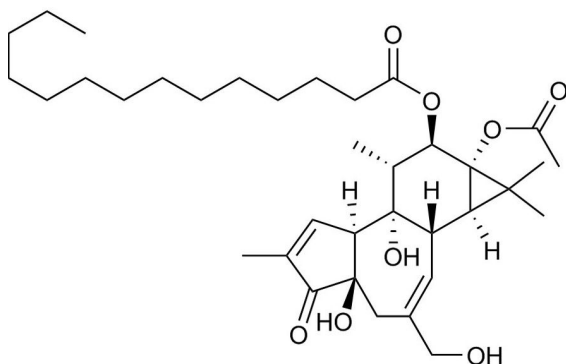

## Properties

|                      |                                                                                                                                                                                                                                                       |
|----------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Physical Appearance: | A crystalline solid                                                                                                                                                                                                                                   |
| Storage:             | Product stable at -20°C as supplied. Protect product from prolonged exposure to light. For long-term storage store with a desiccant. For product expiry date, please contact <a href="mailto:techsupport@stemcell.com">techsupport@stemcell.com</a> . |
| Solubility:          | · DMSO ≤ 40 mM<br>· Absolute ethanol ≤ 40 mM<br>For example, to prepare a 10 mM stock solution in DMSO, resuspend 1 mg in 162 µL of DMSO.                                                                                                             |

Prepare stock solution fresh before use. Information regarding stability of small molecules in solution has rarely been reported, however, as a general guide we recommend storage in DMSO at -20°C. Aliquot into working volumes to avoid repeated freeze-thaw cycles. The effect of storage of stock solution on compound performance should be tested for each application.

Compound has low solubility in aqueous media. For use as a cell culture supplement, stock solution should be diluted into culture medium immediately before use. Avoid final DMSO concentration above 0.1% due to potential cell toxicity.

## Published Applications

### DIFFERENTIATION

- Stimulates cardiac differentiation of mesenchymal stem cells (Seo et al.).
- Promotes hematopoietic differentiation (Clemens et al.).

### IMMUNOLOGY

- Stimulates differentiation of IL-15 stimulated lamina propria CD4<sup>+</sup> T cells into CXCR5<sup>+</sup>CD4<sup>+</sup> cells when used together with ionomycin (Catalog #73722; Sarra et al.).

### CANCER RESEARCH

- In combination with AS101, shows anti-leukemic activity and stimulates human leukemia cell lines to differentiate into macrophage-like cells (Glesne & Huberman; Hayun et al.).

## References

- Adams PD & Parker PJ. (1992) Activation of mitogen-activated protein (MAP) kinase by a MAP kinase-kinase. *J Biol Chem* 267(19): 13135–7.
- Blumberg PM. (1988) Protein kinase C as the receptor for the phorbol ester tumor promoters: Sixth Rhoads Memorial Award Lecture. *Cancer Res* 48(1): 1–8.
- Clemens MJ et al. (1992) The role of protein kinase C isoenzymes in the regulation of cell proliferation and differentiation. *J Cell Sci* 103(4): 881–7.
- Glesne D & Huberman E. (2006) Smad6 is a protein kinase X phosphorylation substrate and is required for HL-60 cell differentiation. *Oncogene* 25(29): 4086–98.
- Hayun M et al. (2007) Synergistic effect of AS101 and Bryostatin-1 on myeloid leukemia cell differentiation in vitro and in an animal model. *Leukemia* 21(7): 1504–13.
- Hurley JH et al. (2008) Taxonomy and function of C1 protein kinase C homology domains. *Protein Sci* 6(2): 477–80.
- Sarra M et al. (2013) IL-15 positively regulates IL-21 production in celiac disease mucosa. *Mucosal Immunol* 6(2): 244–55.
- Seo H-H et al. (2016) The role of nuclear factor of activated T cells during phorbol myristate acetate-induced cardiac differentiation of mesenchymal stem cells. *Stem Cell Res Ther* 7(1): 90.

## Related Small Molecules

For a complete list of small molecules available from STEMCELL Technologies, please visit our website at [www.stemcell.com/smallmolecules](http://www.stemcell.com/smallmolecules) or contact us at [techsupport@stemcell.com](mailto:techsupport@stemcell.com).

**This product is hazardous. Please refer to the Safety Data Sheet (SDS).**

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.