

Cytokines

Human Recombinant EGF, ACF



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #	78136	100 µg
	78136.1	500 µg
	78136.2	1000 µg

Product Description

Epidermal growth factor (EGF) is characterized by high-affinity binding to various EGF receptors (EGFRs) and the production of mitogenic responses (Carpenter & Cohen). EGF promotes EGFR dimerization, resulting in activation of downstream pathways including PI3K, ERK1/2, JAK/STAT, β -catenin, and calcium signaling. EGF is secreted by the gut-associated salivary and Brunner's glands, is found in a variety of body fluids, and stimulates cell proliferation and differentiation in rodent and neonatal human intestine (Wright et al.). Central nervous system stem cells also proliferate in response to the EGF stimulus (Reynolds & Weiss). This product is animal component-free.

Product Information

Alternative Names:	Beta-urogastrone, HOMG4, Epidermal growth factor, Pro-epidermal growth factor, URG, Urogastrone
Accession Number:	P01133
Amino Acid Sequence:	NSDSECPLSH DGYCLHDGVC MYIEALDKYA CNCVVG YIGE RCQYRDLKWW ELR
Predicted Molecular Mass:	6.2 kDa
Species:	Human
Cross Reactivity:	Mouse, Rat, Monkey
Formulation:	Lyophilized from a sterile-filtered aqueous solution containing 0.1% trifluoroacetic acid
Source:	E. coli


Specifications

Activity:	The specific activity is $\geq 1.0 \times 10^7$ units/mg ($EC_{50} \leq 0.1$ ng/mL) as determined by a cell proliferation assay using BALB/c 3T3 cells.
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic Limulus amoebocyte lysate (LAL) analysis and is ≤ 1 EU/µg protein.

Preparation and Storage

Storage:	Store at -20°C to -80°C.
Stability:	Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Resuspend the product in sterile water containing 0.1% bovine serum albumin (BSA) to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. Store at 2 - 8°C for up to 1 month or at -20°C to -80°C for up to 3 months. Avoid repeated freeze-thaw cycles. NOTE: If reconstituted product will be used immediately, BSA is not required.

Data

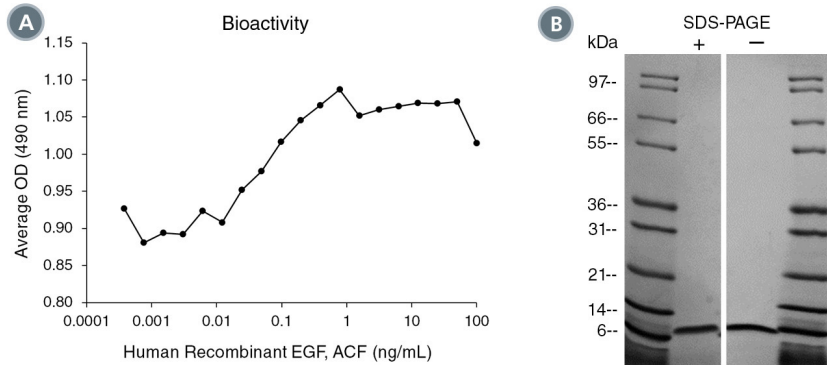

(A) The biological activity of Human Recombinant EGF, ACF was tested by its ability to promote the proliferation of BALB/c 3T3 cells. Cell proliferation was measured using a fluorometric assay method. The EC₅₀ is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC₅₀ in the example above is 0.0459 ng/mL.

(B) 1 µg of Human Recombinant EGF, ACF was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining.

Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Carpenter G & Cohen S. (1990) Epidermal growth factor. *J Biol Chem* 265(14): 7709–12.
- Reynolds BA & Weiss S. (1996) Clonal and population analyses demonstrate that an EGF-responsive mammalian embryonic CNS precursor is a stem cell. *Dev Biol* 175(1): 1–13.
- Wright NA et al. (1990) Induction of a novel epidermal growth factor-secreting cell lineage by mucosal ulceration in human gastrointestinal stem cells. *Nature* 343(6253): 82–5.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.