

Cytokines

Human Recombinant IL-6, ACF



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Interleukin 6, animal component-free

Catalog #	78148	20 µg
	78148.1	100 µg
	78148.2	1000 µg

Product Description

Interleukin 6 (IL-6) is a pleiotropic growth factor with the wide range of biological activities in immune regulation, hematopoiesis, and oncogenesis. IL-6 is produced by a variety of cell types including T cells, B cells, monocytes and macrophages, fibroblasts, hepatocytes, vascular endothelial cells, and various tumor cell lines. On its own or in combination with other factors such as IL-2 and interferon- γ , IL-6 stimulates the proliferation of B cells, T cells, and hybridoma cells (Hirano et al.; Mihara et al.; Tanaka et al.). In combination with cytokines such as IL-3, GM-CSF and SCF, IL-6 has been shown to promote hematopoietic progenitor cell proliferation and differentiation in vitro. IL-6 signals through a cell surface type I cytokine receptor complex consisting of the ligand-binding IL-6 α (CD126) and the signal-transducing gp130 subunits. The binding of IL-6 to its receptor system includes activation of JAK/STAT signaling pathway (Mihara et al.; Peters et al.; Tanaka et al.). This product is animal component-free.

Product Information

Alternative Names:	B cell differentiation factor, BSF-2, IFN- β 2, Interleukin-6
Accession Number:	P05231
Amino Acid Sequence:	MPVPPGEDSK DVAAPHRQPL TSSERIDKQI RYILDGISAL RKETCNKSNM CESSKEALAE NNLNLPKMAE KDGCFQSGFN EETCLVKIIT GLLEFEVYLE YLQNRFESSE EQARAVQMST KVLIQFLQKK AKNLDAITTP DPTTNASLLT KLQAQNQWLQ DMTTHLILRS FKEFLQSSLR ALRQM
Predicted Molecular Mass:	21 kDa
Species:	Human
Cross Reactivity:	Mouse, Rat, Monkey
Formulation:	Lyophilized from a sterile-filtered aqueous solution containing 0.1% trifluoroacetic acid
Source:	E. coli


Specifications

Activity:	The specific activity is $\geq 4.0 \times 10^6$ units/mg ($EC_{50} \leq 0.25$ ng/mL) as determined by proliferation of 7TD1 cells.
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic Limulus amebocyte lysate (LAL) analysis and is ≤ 1 EU/ μ g protein.

Preparation and Storage

Storage:	Store at -20°C to -80°C.
Stability:	Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Resuspend the product in 10 mM hydrochloric acid containing 0.1% bovine serum albumin (BSA) to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. Store at 2 - 8°C for up to 1 month or at -20°C to -80°C for up to 3 months. Avoid repeated freeze-thaw cycles. NOTE: If reconstituted product will be used immediately, BSA is not required.

Data

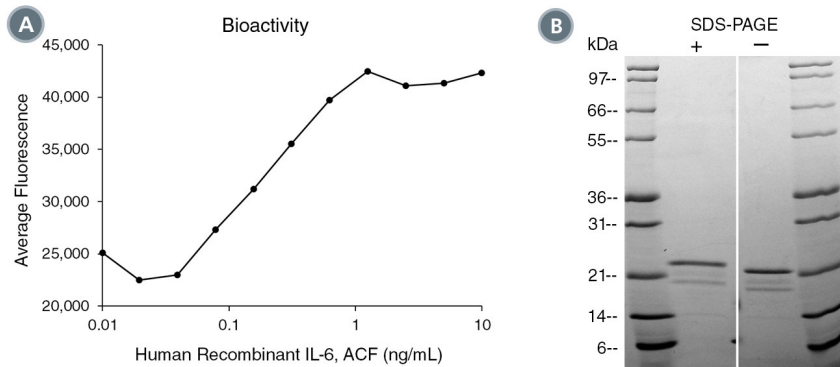

(A) The biological activity of Human Recombinant IL-6, ACF was tested by its ability to promote the proliferation of 7TD1 cells. Cell proliferation was measured after 48 hours of culture using a fluorometric assay method. The EC50 is defined as the effective concentration of the cytokine at which cell proliferation is at 50% of maximum. The EC50 in the example above is 0.2 ng/mL.

(B) 1 µg of Human Recombinant IL-6, ACF was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining.

Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Hirano T et al. Complementary DNA for a novel human interleukin (BSF-2) that induces B lymphocytes to produce immunoglobulin. *Nature* 324(6092): 73–6.
- Mihara M et al. (2012) IL-6/IL-6 receptor system and its role in physiological and pathological conditions. *Clin Sci (Lond)* 122(4): 143–59.
- Peters M et al. (1998) Interleukin-6 and soluble interleukin-6 receptor: direct stimulation of gp130 and hematopoiesis. *Blood* 92(10): 3495–504.
- Tanaka T et al. (2014) IL-6 in inflammation, immunity, and disease. *Cold Spring Harb Perspect Biol* 6(10): a016295.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.