

Cytokines

Human Recombinant PDGF-BB, ACF

Platelet-derived growth factor, animal
component-free



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Catalog #	78154	10 µg
	78154.1	100 µg
	78154.2	1000 µg

Product Description

Platelet-derived growth factor (PDGF) is a dimeric glycoprotein consisting of two disulfide bridge-stabilized polypeptide chains, A and B, which are assembled as heterodimers (PDGF-AB) or homodimers (PDGF-AA and PDGF-BB) (Fretto et al.; Westermark & Heldin). PDGF signals through the receptor tyrosine kinases PDGFRalpha and PDGFRbeta. PDGF-induced migration has been shown to involve MEK/ERK, EGFR, Src, and PI3K/Akt signaling pathways (Kim et al.). PDGF is a potent mitogen for cells of mesenchymal origin- like fibroblasts, glial cells, and vascular smooth muscle cells. PDGF has been implicated in pathogenesis of atherosclerosis, glomerulonephritis, cancer, and in the contraction of vascular smooth muscle cells of rat aortic tissues (Fretto et al.; Sachinidis et al.). PDGF-BB is secreted by osteoblasts to induce mesenchymal stem cell migration and angiogenesis. It has also been shown that PDGF-BB is secreted by preosteoclasts during bone modeling and remodeling to induce angiogenesis and thus proper osteogenesis (Xie et al.). This product is animal component-free.

Product Information

Alternative Names:	GDGF, Glioma-derived growth factor, ODGF, Osteosarcoma-derived growth factor
Accession Number:	P01127
Amino Acid Sequence:	MSLGSLTIAE PAMIAECKTR TEVFEISRRL IDRTNANFLV WPPCVEVQRC SGCCNNRNVQ CRPTQVQLRP VQVRKIEIVR KKPIFKKATV TLEDHLACKC ETVAAARPVT
Predicted Molecular Mass:	12.4 kDa monomer; 24.9 kDa dimer
Species:	Human
Cross Reactivity:	Rat
Formulation:	Lyophilized from a sterile-filtered aqueous solution containing acetic acid
Source:	E. coli


Specifications

Activity:	The specific activity is $\geq 5.0 \times 10^4$ units/mg ($EC_{50} \leq 20$ ng/mL) as determined by proliferation of BALB/c 3T3 cells.
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic Limulus amebocyte lysate (LAL) analysis and is ≤ 1 EU/µg protein.

Preparation and Storage

Storage:	Store at -20°C to -80°C.
Stability:	Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Resuspend the product in sterile water containing 0.1% bovine serum albumin (BSA) to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. Store at 2 - 8°C for up to 1 month or at -20°C to -80°C for up to 3 months. Avoid repeated freeze-thaw cycles. NOTE: If reconstituted product will be used immediately, BSA is not required.

Data

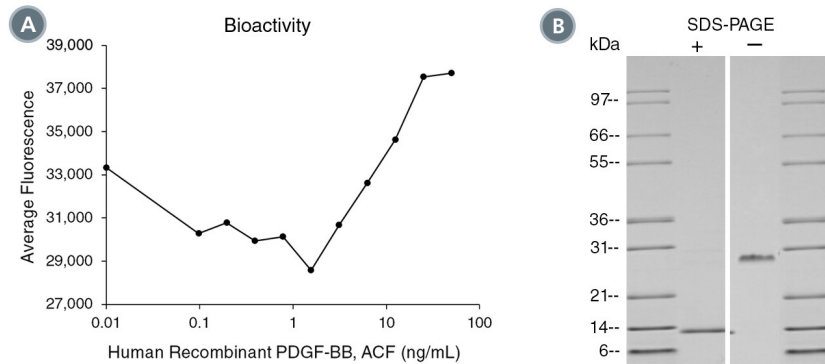

(A) The biological activity of Human Recombinant PDGF-BB, ACF was tested by its ability to promote the proliferation of BALB/c 3T3 cells. Cell proliferation was measured after 24 hours in culture using a fluorometric assay method. The EC₅₀ is defined as the effective concentration of the growth factor at which cell proliferation is at 50% of maximum. The EC₅₀ in the example above is 9.8 ng/mL.

(B) 1 µg of Human Recombinant PDGF-BB, ACF was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining.

Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Fretto LJ et al. (1993) Mechanism of platelet-derived growth factor (PDGF) AA, AB, and BB binding to alpha and beta PDGF receptor. *J Biol Chem* 268(5): 3625–31.
- Kim SJ et al. (2007) Differential effect of FGF and PDGF on cell proliferation and migration in osteoblastic cells. *Growth Factors* 25(2): 77–86.
- Sachinidis A et al. (1990) The platelet-derived growth factor isomers, PDGF-AA, PDGF-AB and PDGF-BB, induce contraction of vascular smooth muscle cells by different intracellular mechanisms. *FEBS Lett* 275(1-2): 95–8.
- Westermarck B & Heldin CH. (1993) Platelet-derived growth factor. Structure, function and implications in normal and malignant cell growth. *Acta Oncol* 32(2): 101–5.
- Xie H et al. (2014) PDGF-BB secreted by preosteoclasts induces angiogenesis during coupling with osteogenesis. *Nat Med* 20(11): 1270–8.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.