

Cytokines

Human Recombinant SCF, ACF



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

| | | |
|-----------|---------|---------|
| Catalog # | 78155 | 10 µg |
| | 78155.1 | 100 µg |
| | 78155.2 | 1000 µg |

Product Description

Stem cell factor (SCF) is an early-acting cytokine that plays a pivotal role in the regulation of embryonic and adult hematopoiesis. SCF promotes cell survival, proliferation, differentiation, adhesion, and functional activation of cells at multiple levels of the hematopoietic hierarchy. Together with other cytokines such as thrombopoietin and Flt3/Flk-2 Ligand, SCF is commonly used to promote expansion of primitive hematopoietic stem cells and multi-potent progenitor cells in culture (Martin et al.; Kent et al.). In synergy with various growth factors, including IL-2, IL-3, IL-6, IL-7, G-CSF, and erythropoietin, SCF increases proliferation and differentiation of myeloid and erythroid progenitor cells and a subset of lymphoid progenitor cells (Broudy). SCF is also a primary growth and activation factor for mast cells and eosinophils. It is produced by stromal cells in the fetal liver, bone marrow, and thymus, in the central nervous system, in keratinocytes, and in the gut mucosa, and can function as a chemotactic and chemokinetic factor.

SCF exists in two biologically active splice forms: a soluble and a transmembrane isoform. Upon binding to its receptor (c-Kit tyrosine kinase receptor; CD117), it activates PI3K, JAK/STAT, and MAPK pathways. SCF and signaling from c-Kit have also been reported to play an important role in pigmentation, fertility, vasculogenesis, motility of the gut via c-Kit positive interstitial cells of Cajal, and in the migration of neuronal stem and progenitor cells to sites of injury in the brain. Signals of proliferation and survival transmitted through c-Kit are also likely to contribute to the initiation and progression of many human malignancies (Lennartsson & Ronnstrand). This product is animal component-free.

Product Information

| | |
|---------------------------|---|
| Alternative Names: | Kit ligand, Mast cell growth factor, Steel factor |
| Accession Number: | P21583 |
| Amino Acid Sequence: | MEGICRNRVT NNVKDVTKLV ANLPKDYMIT LKYVPGMDVL PSHCWISEMV VQLSDSLTDL LDKFSNISEG LSNYSIIDKL VNIVDDLVEC VKENSSKDLK KSFKSPEPRL FTPEEFFRIF NRSIDAFKDF VVASETSDCV VSSTLSPEKD SRVSVTKPFM LPPVA |
| Predicted Molecular Mass: | 18.6 kDa |
| Species: | Human |
| Cross Reactivity: | Weakly active on mouse cells |
| Formulation: | Lyophilized from a sterile-filtered aqueous solution containing acetic acid |
| Source: | E. coli |


Specifications

| | |
|------------------|---|
| Activity: | The specific activity is $\geq 6.7 \times 10^4$ units/mg ($EC_{50} \leq 15$ ng/mL) as determined by a cell proliferation assay using TF-1 cells. |
| Purity: | $\geq 95\%$ |
| Endotoxin Level: | Measured by kinetic Limulus amoebocyte lysate (LAL) analysis and is ≤ 1 EU/µg protein. |

Preparation and Storage

| | |
|--------------|---|
| Storage: | Store at -20°C to -80°C. |
| Stability: | Stable as supplied for 12 months from date of receipt. |
| Preparation: | Centrifuge vial before opening. Resuspend the product in 10 mM hydrochloric acid containing 0.1% bovine serum albumin (BSA) to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex. Store at 2 - 8°C for up to 1 month or at -20°C to -80°C for up to 3 months. Avoid repeated freeze-thaw cycles. NOTE: If reconstituted product will be used immediately, BSA is not required. |

Data

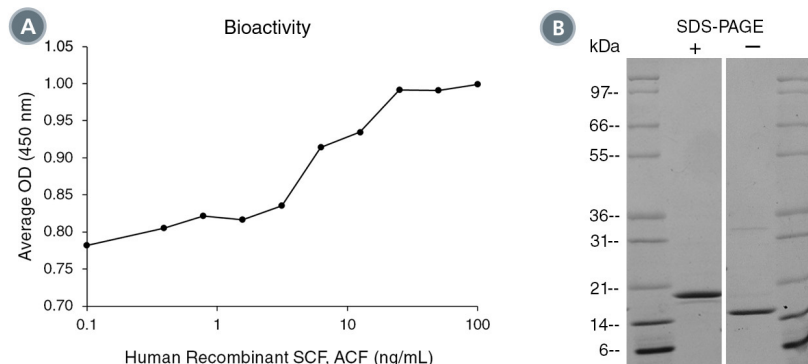

(A) The biological activity of Human Recombinant SCF, ACF was tested by its ability to promote the proliferation of TF-1 cells. Cell proliferation was measured after 72 hours of culture using a fluorometric assay method. The EC50 is defined as the effective concentration of the cytokine at which cell proliferation is at 50% of maximum. The EC50 in the example above is 5.3 ng/mL.

(B) 1 µg of Human Recombinant SCF, ACF was resolved with SDS-PAGE under reducing (+) and non-reducing (-) conditions and visualized by Coomassie Blue staining.

Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, please visit our website at www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Broudy VC. (1997) Stem cell factor and hematopoiesis. *Blood* 90(4): 1345–64.
- Kent D et al. (2008) Regulation of hematopoietic stem cells by the steel factor/KIT signaling pathway. *Clin Cancer Res* 14(7): 1926–30.
- Lennartsson J & Ronnstrand L. (2012) Stem cell factor receptor/c-Kit: from basic science to clinical implications. *Physiol Rev* 92(4): 1619–49.
- Martin FH et al. (1990) Primary structure and functional expression of rat and human stem cell factor DNAs. *Cell* 63(1): 203–11.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2017 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.