

Cytokines

Human Recombinant IL-31, ACF



Scientists Helping Scientists™ | WWW.STEMCELL.COM

TOLL FREE PHONE 1 800 667 0322 • PHONE +1 604 877 0713

INFO@STEMCELL.COM • TECHSUPPORT@STEMCELL.COM

FOR GLOBAL CONTACT DETAILS VISIT OUR WEBSITE

Interleukin 31, animal component-free

Catalog #	78217	25 µg
	78217.1	100 µg
	78217.2	500 µg
	78217.3	1000 µg

Product Description

Interleukin 31 (IL-31), a four-helix bundle inflammatory cytokine, belongs to the IL-6 cytokine family which includes IL-6 (Catalog #78050), oncostatin M (Catalog #78094), leukemia inhibitory factor (LIF; Catalog #78055), and cardiotrophin-1 (Catalog #78160; Dillon et al.). IL-31 signals through a heterodimer composed of IL-31RA (also known as gp130-like receptor or GPL) and the oncostatin-M receptor (OSMR), both of which are expressed on monocytes (Diveu et al.; Ghilardi et al.), epithelial cells (Ip et al.), and keratinocytes (Kato et al.). Signaling through the GPL/OSMR complex activates the JAK/STAT, RAS/ERK, and PI3K/AKT signaling pathways, resulting in the downstream activation of STAT1, STAT3, and STAT5 transcription factors (Cornelissen et al.; Dambacher et al.; Dillon et al.; Ip et al.). IL-31 responses have been associated with allergic responses and inflammatory skin diseases including atopic dermatitis (Cornelissen et al.; Gangemi et al.). This product is animal component-free.

Product Information

Alternative Names:	Interleukin-31
Accession Number:	Q6EBC2
Amino Acid Sequence:	SHTLPVRLLR PSDDVQKIVE ELQSLSKMLL KDVEEEKGVL VSQNYTLPC L SPDAQPPNNI HSPAIRAYLK TIRQLDNKSV IDEIIEHLDK LIFQDAPETN ISVPTDTHC KRFILTISQQ FSECMDLALK SLTSGAQQAT T
Predicted Molecular Mass:	15.8 kDa
Species:	Human
Cross Reactivity:	Not determined
Formulation:	Lyophilized from a sterile-filtered solution containing sodium phosphate and sodium chloride, pH 7.5.
Source:	E. coli

Specifications


Activity:	Biological activity was detectable at ≤ 25 ng/mL as determined by a STAT3 activation assay using U-87 MG cells.
Purity:	$\geq 95\%$
Endotoxin Level:	Measured by kinetic Limulus amebocyte lysate (LAL) analysis and is ≤ 1 EU/µg protein.

Preparation and Storage

Storage:	Store at -20°C to -80°C.
Stability:	Stable as supplied for 12 months from date of receipt.
Preparation:	Centrifuge vial before opening. Reconstitute the product in sterile water to at least 0.1 mg/mL by pipetting the solution down the sides of the vial. Do not vortex.

OPTIONAL: After reconstitution, if product will not be used immediately, dilute with concentrated bovine serum albumin (BSA) to a final BSA concentration of 0.1%. The effect of storage of stock solution on product performance should be tested for each application. As a general guide, do not store at 2 - 8°C for more than 1 month or at -20°C to -80°C for more than 3 months. Avoid repeated freeze-thaw cycles.

Data

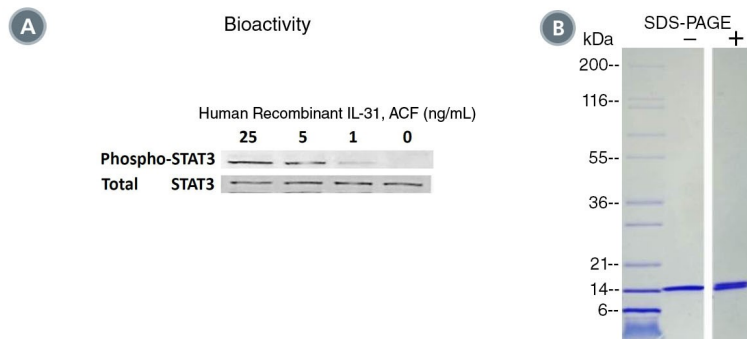

(A) The biological activity of Human Recombinant IL-31, ACF was tested by STAT3 activation in U-87 MG cells. The detection limit was ≤ 25 ng/mL.

(B) Human Recombinant IL-31, ACF was resolved with SDS-PAGE under reducing (+) conditions and non-reducing (-) conditions and visualized by Coomassie Blue staining. Human Recombinant IL-31, ACF has a predicted molecular mass of 15.8 kDa.

Related Products

For a complete list of cytokines, as well as related products available from STEMCELL Technologies, visit www.stemcell.com/cytokines or contact us at techsupport@stemcell.com.

References

- Cornelissen C et al. (2012) Signaling by IL-31 and functional consequences. *Eur J Cell Biol* 91(6–7): 552–66.
- Dambacher J et al. (2007) Interleukin 31 mediates MAP kinase and STAT1/3 activation in intestinal epithelial cells and its expression is upregulated in inflammatory bowel disease. *Gut* 56(9): 1257–65.
- Dillon SR et al. (2004) Interleukin 31, a cytokine produced by activated T cells, induces dermatitis in mice. *Nat Immunol* 5(7): 752–60.
- Diveu C et al. (2003) GPL, a novel cytokine receptor related to GP130 and leukemia inhibitory factor receptor. *J Biol Chem* 278(50): 49850–9.
- Gangemi S et al. (2017) Interleukin 31 and skin diseases: A systematic review. *Allergy Asthma Proc* 38(6): 401–8.
- Ghilardi N et al. (2002) A novel type I cytokine receptor is expressed on monocytes, signals proliferation, and activates STAT-3 and STAT-5. *J Biol Chem* 277(19): 16831–6.
- Ip WK et al. (2007) Interleukin-31 induces cytokine and chemokine production from human bronchial epithelial cells through activation of mitogen-activated protein kinase signalling pathways: implications for the allergic response. *Immunology* 122(4): 532–41.
- Kato A et al. (2014) Distribution of IL-31 and its receptor expressing cells in skin of atopic dermatitis. *J Dermatol Sci* 74(3): 229–35.

STEMCELL TECHNOLOGIES INC.'S QUALITY MANAGEMENT SYSTEM IS CERTIFIED TO ISO 13485. PRODUCTS ARE FOR RESEARCH USE ONLY AND NOT INTENDED FOR HUMAN OR ANIMAL DIAGNOSTIC OR THERAPEUTIC USES UNLESS OTHERWISE STATED.

Copyright © 2019 by STEMCELL Technologies Inc. All rights reserved including graphics and images. STEMCELL Technologies & Design, STEMCELL Shield Design, and Scientists Helping Scientists are trademarks of STEMCELL Technologies Canada Inc. All other trademarks are the property of their respective holders. While STEMCELL has made all reasonable efforts to ensure that the information provided by STEMCELL and its suppliers is correct, it makes no warranties or representations as to the accuracy or completeness of such information.